

G3D Animal Monitor Operating Manual

Version V1.0

GENERAL MEDITECH, INC.

December 24, 2007

Information in this document is subject to change without notice. No part in this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of GENERAL MEDITECH, INC.

Contact:

South Office 4/F, Kezhi Rd.W.No. 1, Science Park,
Nanshan District, Shenzhen, Guangdong P.R.China

Tel: 86-755-26500832

Fax: 86-755-26546285

Zip: 518057

[Http://www.szmedtech.com](http://www.szmedtech.com)

E-mail:zhonglin@szonline.net

CONTENTS

CHAPTER 1 OVERVIEW	1
1.1 Brief introduction.....	1
1.2 Warranty.....	1
1.3 The equipment installation.....	1
1.3.1 Opened package inspection.....	1
1.3.2 Connect AC power and network	2
1.3.3 Other requirements on the connection of external equipments	2
1.4 Safety information	2
1.4.1 The electrical safety classifications.....	2
1.4.2 Prohibition.....	2
1.4.3 Cautions	3
1.5 Function composition	3
1.6 Product indexes.....	3
1.6.1 Performance index	3
1.6.2 The available environment conditions for the monitor	5
CHAPTER 2 INTRODUCTION	6
2.1 External appearance.....	6
2.1.1 Front panel	6
2.1.2 Rear panel.....	8
2.2 Introduction on display	9
CHAPTER 3 BASIC OPERATION	11
3.1 How to use rotating mouse	11
3.2 How to use the volume knob	11
3.3 How to use keys.....	11
3.4 Menu operation.....	11
3.5 Set date and time.....	12
3.6 Setting unit.....	12
3.7 Change the waveform speed.....	12
3.8 Setting monitoring format.....	12
3.9 Setting Alarm.....	12
3.10 Preparatory work	15
3.11 Remote maintenance.....	15
3.11.1 Setting your monitor.....	15
CHAPTER 4 ECG MONITORING	25
4.1 ECG electrode monitoring.....	25
4.1.1 Installing the ECG electrode	26
4.1.2 ECG monitoring interface	27
4.2 Parameter display.....	29
4.3 Maintenance.....	30
4.4 Trouble shooting	30
CHAPTER 5 NIBP MONITORING.....	31

5.1 Overview	31
5.2 Blood pressure cuff installation	31
5.3 The limitations on the measurement	32
5.4 Normal Blood pressure values	33
5.5 Monitoring display interface.....	34
5.6 Introduction on functions.....	34
5.6.1 Blood pressure measurement	34
5.6.2 Venipuncture	35
5.7 Maintenance and cleaning	36
5.8 Trouble shooting	37
CHAPTER 6 THE IMPEDANCE RESPIRATION MONITORING.....	38
6.1 The electrode connection	38
6.2 RESP monitoring inspection.....	38
6.3 Monitoring interface	39
6.4 Prompting messages	40
CHAPTER 7 NASAL TUBE RESPIRATION MONITORING.....	41
7.1 Monitoring interface	41
7.2 The status prompting messages	41
CHAPTER 8 TEMPERATURE MONITORING	42
8.1 The installation of the temperature sensor	42
8.2 Body temperature parameter menu	43
8.3 Maintenance and cleaning	43
CHAPTER 9 IBP SETTING AND MEASURING (OPTIONAL)	44
9.1 Waveforms displaying	44
9.2 Module status display	45
9.3 Pressure value display.....	45
CHAPTER 10 TREND ANALYSIS.....	48
10.1 Main interface trend analysis.....	48
10.2 The trend analysis of the trend window	50
CHAPTER 11 THE WAVEFORM FREEZE AND REPLAY	53
11.1 Waveform freeze operational specification.....	53
CHAPTER 12 ENVENT RECORD.....	54
CHAPTER 13 HOW TO USE THE DRUG CALCULATOR.....	55
CHAPTER 14 THE PRINTING FUNCTION OF THE MONITOR (OPTIONAL).....	61
14.1 The connection mode of the printer	61
14.2 The configuration of the monitor	61
14.3 Waveform and parameter printing	62
CHAPTER 15 CO2 SETTING AND MEASURING (OPTIONAL).....	63
15.1 Introduction	63

15.2 CO2 measuring	63
15.3 Monitoring interface	64
CHAPTER 16 PACKAGING AND LABEL	65
16.1 Packaging frontispiece.....	65
16.2 Packaging profile.....	65
16.3 Product nameplate	65

Chapter 1 Overview

Welcome to use G3D Multi-Parameter Animal Monitor! The main purpose of this Manual is to provide operating guide, instrument repair and maintenance information to users. The detailed performance indexes, routine installation, operation and maintenance method as well as safety information are illustrated in this Manual. Before use, the operator shall carefully read this Manual in order to properly and correctly operate this monitor, so that it can reach specified safety standards and performance indexes.

1.1 Brief introduction

The G3D Multi-Parameter Animal Monitor is a multifunctional overall physiologic index monitor, it is applicable to monitor the ECG, NIBP, SpO₂, HR, RESP, TEMP in hospitals.

This product has the multi-parameter functions which can be selected, combined and configured by the user according to his requirements (but the G3D Multi-Parameter Animal Monitor which you are using probably only had part parameters and corresponding accessories selected by you when you bought it).

This product is composed of main machine and corresponding functional accessories such as ECG cable, blood pressure cuff, SpO₂ sensor, body temperature sensor and BP hose.

This product has four input/output interfaces for printer, network communication, external CRT and external telephone line.

1.2 Warranty

Our Company offers a 1.5 year warranty of free charge for the G3D Multi-Parameter Animal Monitor since the date of the purchase, while offers 1 year warranty for the SpO₂ sensor, blood pressure cuff, ECG cable and body temperature sensor since the date of the purchase, the consumables are excluded from the scope of this warranty.

The following situations are not included in the scope of the warranty:

The serial number of the monitor was ripped off or unreadable.

The monitor was damaged due to improper connection with other equipment.

The monitor was damaged due to accident.

The user modified the monitor without our Company's written authorization.

1.3 The equipment installation

1.3.1 Opened package inspection

Please open the package according to the positions marked on the package box, and take out the monitor and accessories carefully.


- Count the accessories according to the packing list,
- Check the monitor and accessories if any damages occurred, e.g. mechanical deterioration.

If has any problem, please contact the distributor immediately.

Friendly remind: The package box and protective filling material shall be well preserved for your further transportation and storage.

1.3.2 Connect AC power and network

- Confirm the AC power source is in conformity with the requirements of this equipment: (90~256)V a.c., 50/60Hz;
- Connect the monitor to the AC power source with the flexible cable adhered to the monitor.
- The connection with the external equipments shall be in conformity with the requirements specified in the Manual.
- Electric potential balance terminal (refer to the rear panel introduction) shall be connected to the protective grounding end of the public power grid or properly connected to the other protective grounding ends.

 If you found and proved that the functions of the monitor have been damaged or it shown error prompting message, the monitor can not be used for monitoring on the animal, and you shall contact the relevant department in hospital or manufacturer.

1.3.3 Other requirements on the connection of external equipments



The parallel printer faucet, network faucet are on the rear panel of the monitor, through which the monitor can be connected to the external printer and can be connected to the network to update the monitor to add more new functions.

When the monitor is connected to the external equipment via the faucets of this equipment, please refer to GB 9706.15.


Do not place the monitor against the wall, and do not plug up the heat emission holes on the back and two sides of the monitor.

1.4 Safety information

1.4.1 The electrical safety classifications


The electrical safety classification of this product is I class internal stand-by power equipment; if there is a symbol  (electrocardiogram) beside the sensor insertion slot, its application section is CF defibrillation prevention type, if there is a symbol  beside the sensor insertion slot, its application section (except the ECG application section) is BF defibrillation prevention type.

The marked symbol  on the Manual means: attention, refer to document attached with the machine.

Anywhere marked with , means when you using this equipment, you shall pay more attention on it in order to guarantee the safety of the animal, operator or equipment.

1.4.2 Prohibition

 Do not use the G3D Multi-Parameter Animal Monitor for asphyxiation monitoring.

 Do not use the G3D Multi-Parameter Animal Monitor during the magnetic resonance imaging (MRI) or CT inspection process.

⚠ Do not use this monitor in the environment filled with inflammable gas and anaesthesia gas.

1.4.3 Cautions

⚠ In order to make the monitor safe grounding, the hospital must provide the power socket with complete wire under voltage, zero line and protective grounding in accordance with national standard. Or the hospital bears all consequences.

⚠ Do not dismantle the housing of the monitor at will for the sake of electric shock. The maintenance and upgrade of the monitor shall be carried out by the qualified and authorized service personnel.

⚠ Only qualified accessories can be used for the monitor.

⚠ Keep the monitor dry, prevent it from water and humidity, and avoid strong collision.

⚠ In order to make the monitor safe grounding, hospital must provide the power socket with complete wire under voltage, zero line and protective grounding.

⚠ Do not use this monitor in the environment filled with inflammable gas and anaesthesia gas.

⚠ Do not dismantle the monitor cover board for the sake of electric shock.

⚠ When using the electrosurgery unit during the monitoring, the loop of the electrosurgery unit shall be properly connected to prevent from burning, even the life risk.

⚠ When using this monitor, the F-type application section cannot be connected to other conductivity or ground.

⚠ In order to avoid losing time for diagnosis and treatment, please configure adequate alarm settings according to different conditions of each animal (for detailed setting methods, refer to the 3.9 Setting Alarm and 3.11 Alarm menu in this Manual.

⚠ When various equipments are connected to the same animal, pay attention to the danger of leakage current overlap.

⚠ Do not put the monitor in environments such as high temperature and high pressure, gas fumigation or liquid immersion. Before clean or sterilize the monitor, cut off the general power please.

⚠ The ECG Cable uses five lead cable only, and cannot be connected to other signal terminals.

⚠ G3D Multi-Parameter Animal Monitor has no defibrillation synchronization, so it cannot be connected to the synchronization defibrillation instrument.

1.5 Function composition

The G3D Multi-Parameter Animal Monitor can be used to monitor the major parameters such as electrocardiogram (ECG), body temperature (TEMP), respiration (RESP), saturation oxygen (SpO₂), invasive blood pressure (IBP), noninvasive blood pressure (NIBP). It integrates the parameter measurement module function, display and record export as a whole, which forms a compact and portable monitor.

1.6 Product indexes

1.6.1 Performance index

◆ Electrocardiogram (ECG)

Heart rate (HR)

Measurement range: 30BPM ~ 350BPM;

Accuracy in measurement: 1BPM or $\pm 1\%$, the larger one is applicable.

Heart rate alarm settings and allowable deviation:

upper limit (60 ~ 300)bpm, allowable deviation $\pm 10\%$

lower limit (10 ~ 200)bpm, allowable deviation $\pm 10\%$

ECG noise level: $\leq 30\mu\text{V}$

ECG frequency response: filter mode: (1 ~ 25)Hz

Non-filter mode: (0.05 ~ 100)Hz

ECG input loop current: $\leq 0.1\mu\text{A}$

ECG display sensitivity: continuous adjustment

Scan speed grades and errors:

The scan rate can be divided into four grades for option: 6.25mm/s, 12.5mm/s, 25mm/s, and 50mm/s.

The scan rate error: $\leq \pm 10\%$

◆ Noninvasive blood pressure (NIBP)

Measurement range

a) Systolic blood pressure (SBP): (5.3 ~ 36.0) kPa, or (40 ~ 270) mmHg

b) Diastolic blood pressure (DBP): (2.7 ~ 26.7) kPa, or (20 ~ 200) mmHg

c) Mean pressure (MAP): (4.0 ~ 29.3) kPa, or (30 ~ 220) mmHg

Accuracy in measurement: $\leq \pm 0.4$ kPa (3 mmHg) or $\pm 2\%$, the larger one is applicable.

◆ Respiration rate (RR)

Measurement range:

0 BPM ~ 120 BPM;

Accuracy in measurement:

± 1 BPM or $\pm 5\%$, the larger one is applicable.

◆ Body temperature (TEMP)

Body temperature channel

The number of channel: 2.

Body temperature sensor configuration: one intra-cavity probe, one optional body surface sensor.

Measurement range: 25 °C ~ 45 °C

Accuracy in measurement: ± 0.2 °C

Response time: ≤ 3 min

◆ Pulse saturation oxygen (SpO₂)

Measurement range: (0 ~ 100) %

Accuracy in measurement: within the range of (70 ~ 100) % ± 2 %

Within the range of (50 ~ 69) % ± 3 %

◆ Pulse rate (PR)

Measurement range:

(30 ~ 250) bpm;

Accuracy in measurement:

1BPM or $\pm 2\%$, the larger one is applicable.

◆ Trend retrospection: it can retrospect the data of past 360 hours.

1.6.2 The available environment conditions for the monitor

Temperature range: working (+5 ~ +45)°C
Transportation and storage (-20 ~ +55)°C

Relative humidity: working ≤ 80 %
Transportation and storage ≤ 93 %

Air pressure: working (700 ~ 1060)hPa
Transportation and storage (500 ~ 1060)hPa

Adaptive power: (90~256)V a.c. 50/60Hz
Maximum input power 80VA
Fuse: T 1.6 A

Continuous working time:

Power working time: ≥ 8 hours

Internal power continuous working time: > 1 hour, (when the battery is fully charged and it measures the blood pressure every ten minutes).

Signal Faucet:

CRT (external monitor):	Standard VGA faucet (optional)
PRINTER (external printer)	PC/AT parallel faucet (optional)
NET (external central monitor)	RJ45 faucet
LINE (external telephone line)	MODEM faucet

Chapter 2 Introduction

2.1 External appearance

2.1.1 Front panel

The front panel of the G3D Multi-Parameter Animal Monitor is shown as Figure 2.1:

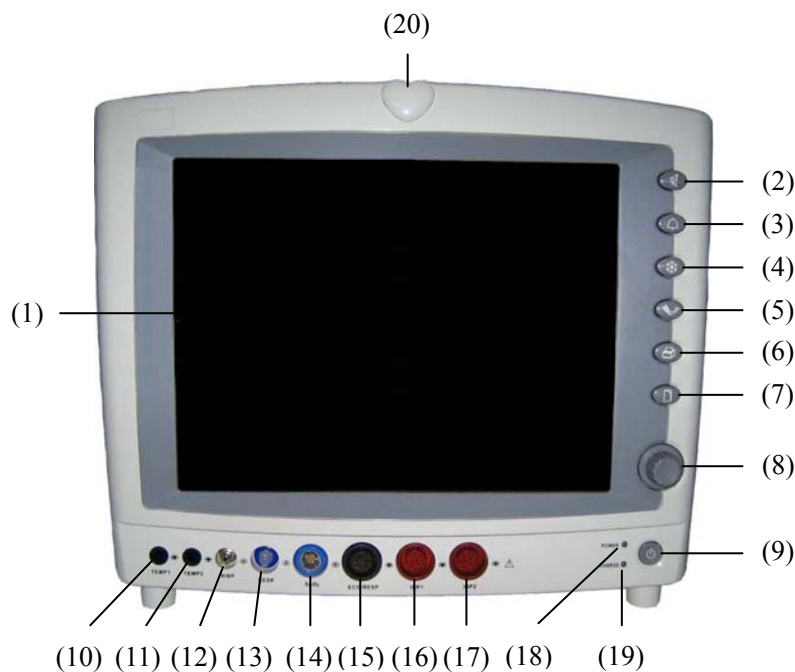






Figure 2.1


(1) **DISPLAY**: Waveforms, menu, alarm and measuring parameters are displayed here.


(2)  — Silence key. Press this button to start or stop the monitor tone.

(3)  — Alarm key. Press this button to enter alarm setting.

(4)  — Freeze key. When the waveforms are updated, press this button to freeze the displayed waveforms.

(5)  — NIBP key. In hold condition, press this button to start NIBP measurement. If the monitor is performing a NIBP measurement, press this button to stop current NIBP measurement.

(6)  — Print key. When the device is equipped with a recording/printer-device, press this button to activate printout of the monitoring results. When the printer is operating, stop print by pressing this button.

(7)  — Main menu key. When no menu is shown in display, press this button to show the system's main menu. Whatever status the menu appears in, exit the menu by pressing this button.

(8) **ROTATING MOUSE**: Rotate the mouse clockwise or counterclockwise to select the menu item. Press the mouse to operate the menu item.

(9) **Power switch**: Press this key for more than two seconds to turn ON and OFF the monitor.

- (10) TEMP1 — Temperature probe connector (channel 1)
- (11) TEMP2 — Temperature probe connector (channel 2)
- (12) NIBP — NIBP cuff hose connector
- (13) RESP — Respiration pipe connector
- (14) SpO₂ — SpO₂ probe connector
- (15) ECG/RESP — ECG cable connector
- (16) IBP1 — IBP transducer connector (channel 1)
- (17) IBP2 — IBP transducer connector (channel 2)
- (18) POWER — Power indicator

ON: The monitor power is turned on.


OFF: The monitor power is turned off.

- (19) CHARGE — Charge indicator

ON: AC power is applied to the monitor.

OFF: AC power is not applied to the monitor.

- (20) Alarm Indicator. When a parameter value exceeds its alarm limit, this indicator flashes once a second.

 Type CF applied part. The unit displaying this symbol contains an F-type isolated (floating) animal part providing a high degree of protection against shock, and is suitable for use during defibrillation.

 TYPE BF applied part. Defibrillator-proof protection against electrical shock.

 Attention: Consult accompanying documents (this manual).

2.1.2 Rear panel

The rear panel of the G3D Multi-Parameter Animal Monitor is shown as Figure 2.2:

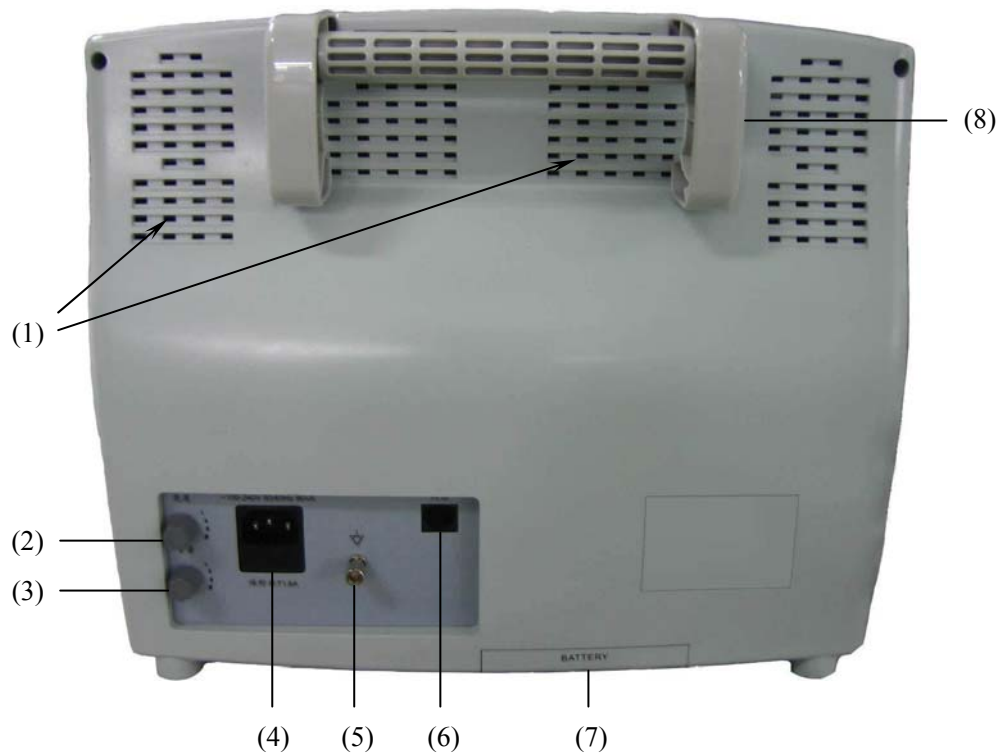


Figure 2.2

(1) Dispersion holes

(2) Brightness adjustment knob

Turn this knob to adjust the brightness of the monitor screen. Turn clockwise to brighter the screen; turn counterclockwise to darker the screen.

(3) Volume adjustment knob

Turn this knob to adjust the volume of the monitor. Turn clockwise to louder the volume; turn counterclockwise to lower the volume.

(4) AC power input connector

A three-wire power cord can be connected to this receptacle to provide AC power supply to the animal monitor.

(5) Electric potential balance terminal

It is marked with symbol ∇ , when the monitor is used with other equipment, it can be connected to equal the voltage.

(6) Network connector: Standard RJ45 connector



Through network connector, this monitor can be connected with the central monitoring system, another monitor, or a PC. It enables the functions of other animal viewing, data output and on-line software upgrading.

(7) Battery cover

(8) Handle

2.2 Display

The display of this monitor is divided into four major areas such as channel display area, parameter display area, status bar and menu bar. The status bar is displayed on the top part of the screen; the channel display area is displayed under the status bar, the left part of the screen, while the menu bar is always displayed at the low part of the screen. The parameters are always displayed at the right part of the screen. The menu bar is a popup one, when is popup, it will shut out part of the channel display area and parameter display area. The basic layout of the display is as following:

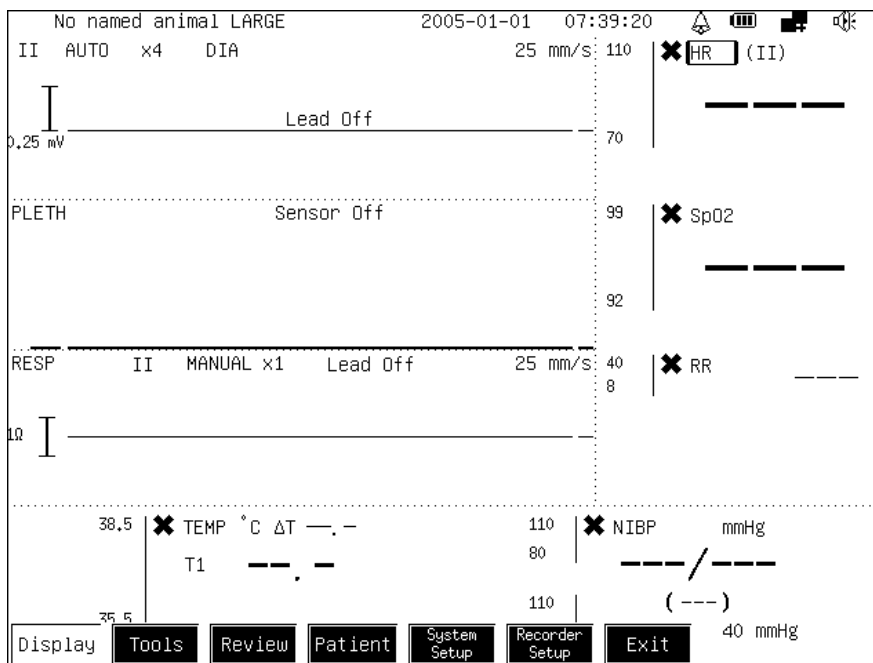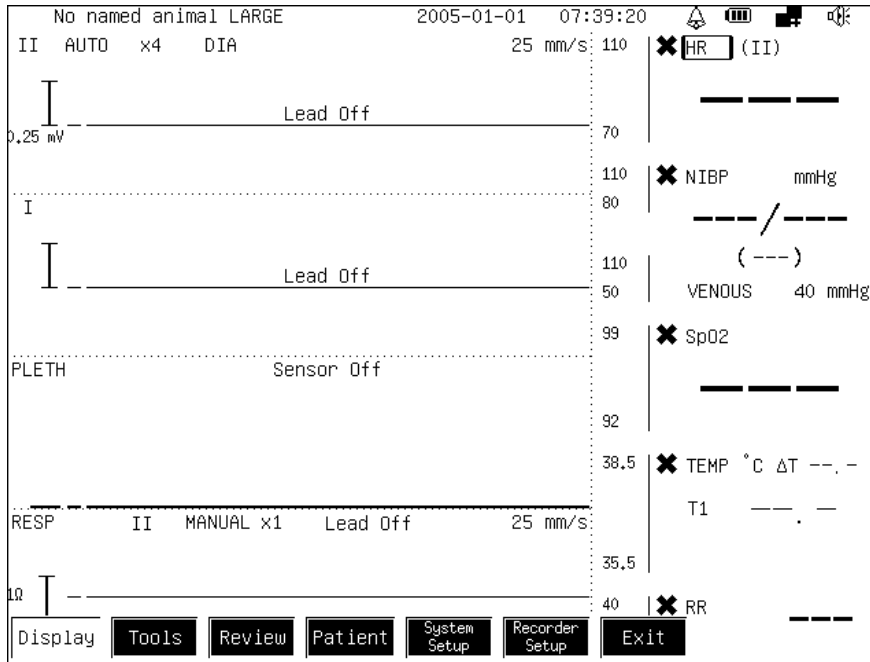

Because the Monitor has provided the channel display area configuration function, the display content at the channel display area can be set by the operator. The parameter display area can be automatically adjusted according to the channel configuration format, so that each parameter can be equally displayed at the right part or right low part of the screen along with corresponding waveform channels.

Channel display area

The channel display area is divided into several channels. Total seven channels are provided; the content of the channel can be selected from any of them such as the ECG waveform, RESP waveform, blood volume waveform, trend graph and trend table. The name of each channel is usually displayed on the left corner of the channel, select the channel, you can visit the menu of the channel. What is displayed at the same row of the channel name is the status information concerning this channel.


Parameter display area







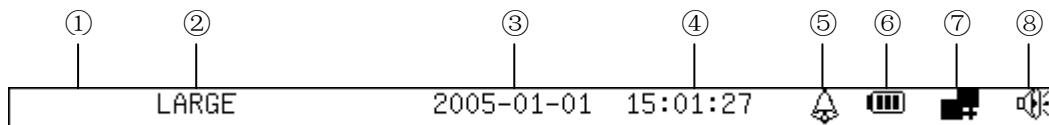

The heart rate (HR), oxygen saturation (SpO2), respiration rate (RESP), body temperature (T1 / T2) and noninvasive blood pressure (NIBP) are displayed in this area, every parameter has a corresponding label; The important parameters are configured with alarm prohibition signs, the prohibition signs are displayed on the left part of the label, the operational menu of relevant user's parameters can be gained by selecting parameter label.


Menu bar

The menu bar area is dynamic popup and disappearance, when no menu is displayed on the screen, the main menu key, the alarm key and the freeze key on the front panel can activate the menus to display, the relevant menus can also be activated by selecting the labels on the display with the rotating mouse.

Status bar



- ① Information bar. Display the Monitor status information and animal's information.
- ② Animal types: large or small.
- ③ Date. Its display format can be modified in the menu.
- ④ Time. Its display format can be modified in the menu.
- ⑤ The alarm status icon:  means the alarm is allowed,  means the alarm is prohibited.
- ⑥ The battery electrical voltage icon: denotes the remaining charge in the battery. Three grids  means full charge; blank  means to be exhausted, charge in time please.
- ⑦ Central monitor network function status icon. Red icon  means the monitor is not connected to the central monitor system, the blue icon  means the monitor has been connected to the central monitor system.
- ⑧ The mute status icon:  means the heartbeat sound is unclosed,  means the heartbeat sound is closed.

⚠ When the symbol  shows, the system cannot produce alarm sound prompting, so the operator shall be especially cautious to use this function.

Chapter 3 Basic Operation

3.1 How to use rotating mouse

At the main interface, there is a yellow frame which means your current position, you can select different sign names with the rotating mouse, or you can click at the current sign name to popup the corresponding functional menu.

3.2 How to use the volume knob

By tuning the volume knob you can adjust the volume, tune counterclockwise, the volume becomes low, tune clockwise, the volume becomes loud.

3.3 How to use keys

- Silence key: turn the heartbeat sound on or off.
- Alarm key: popup or close the alarm menu.
- Freeze key: popup or close the freeze menu.
- NIBP key: inflate the air when measuring the blood pressure.
- Print key: print the waveforms.
- Menu key: popup or close the menu.

3.4 Menu operation

Pop-up the menu:

Press down the corresponding menu keys (menu, freeze, alarm key): the corresponding menu will pop up.

Select the concerned waveform or parameter label on the screen by the rotating mouse, and then click the rotating mouse, the corresponding menu pops out.

Menu browse:

In the menu bar, the knob of black background and white character is the knob selected by the cursor. Using the rotating mouse to select the required knob, then press the rotating mouse, and you can get the function of this knob.


Return to the previous menu:

Select the “Exit” knob at the right end of the menu bar, the menu bar returns to the previous menu.

If press the menu key on the front panel when in the menu displaying status, it will exit all the menus.

If no key operation for more than one minute, the Monitor will automatically exit all the menus.

Main menu:

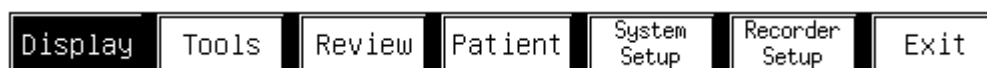

3.5 Set date and time

Note: Changing date or time will influence the storage of the trend and events, after the modification, you must turn off the Monitor and turn it on again immediately!

Select menu key -> System Setup -> Set Time



Set Date


Year: 2005 Month: 1 Day: 1

Hour: 1 Minute: 1 Second: 15

Set system time will disturb trend

3.6 Setting unit

Select menu key -> System Setup -> Region Settings -> Select Units

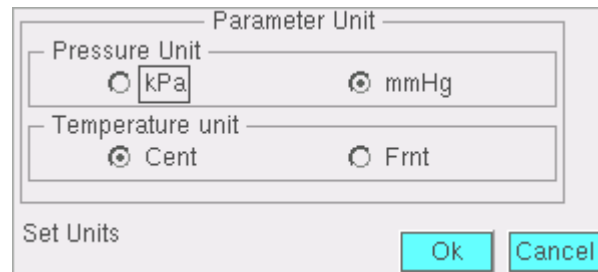

Parameter Unit

Pressure Unit

kPa mmHg

Temperature unit

Cent Frnt

Set Units

3.7 Change the waveform speed

Decrease the waveform speed and the waveform is compressed, you will see longer time quantum, increase the waveform speed, the waveform is expanded, you get more detailed view, so you can change the ECG waveform speed and RESP waveform speed.


3.8 Setting monitoring format

The common monitoring formats include 1 ECG, 2 ECGs, 3 ECGs, 7 ECGs and Huge Digit, and the other five formats defined by the user. The detailed information for the format you required, please refer to the Format Setup.


3.9 Setting Alarm

All the parameters on the screen can be set with the alarm and alarm limit, you can use the rotating mouse to select the parameter you want to set the alarm and then click the mouse to operate the relevant menu.

- Open or close the alarm
 - a. Open or close the alarm of a certain parameter
 - Select the sign name of the parameter > Press the mouse -> Select the alarm on/off menu**

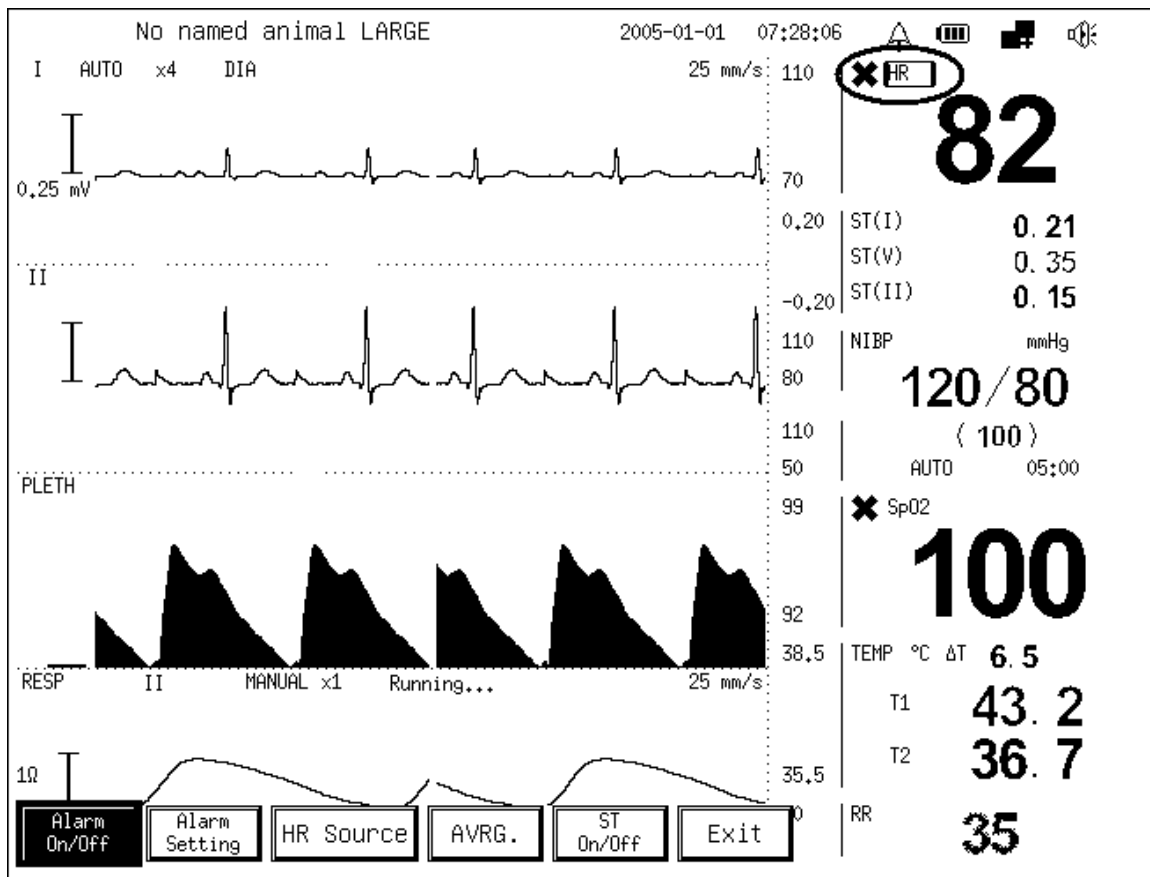


- b. Open or close the alarm of all parameters
Select **alarm key** -> **Alarm On/Off**

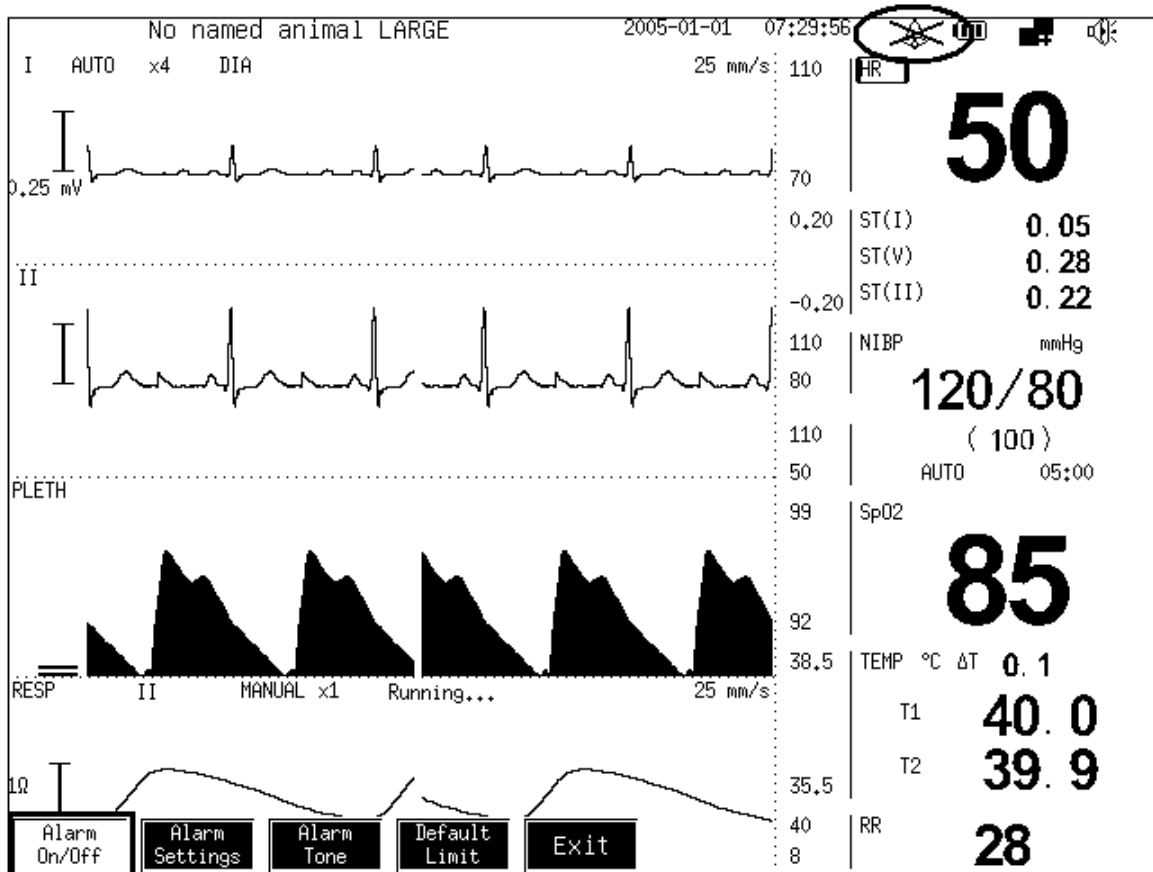


- Alarm setting
 - The alarm setting for a certain parameter

Select the sign name of this parameter -> Press the mouse -> Select the alarm setting menu

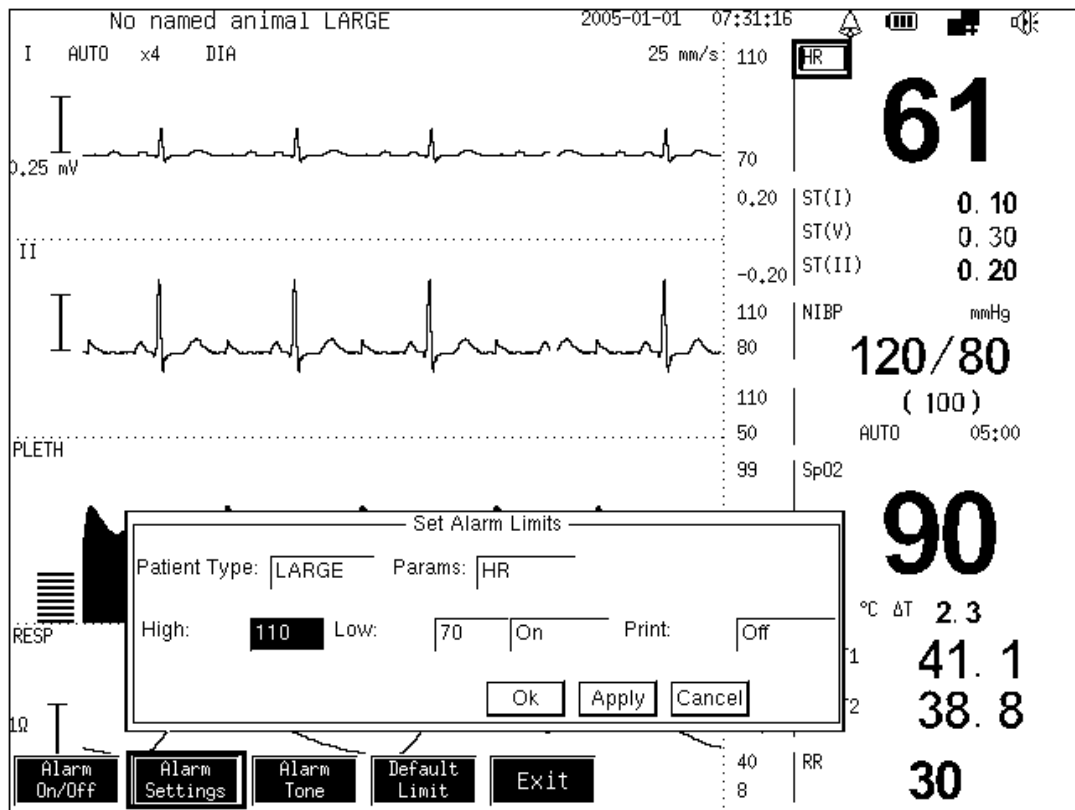

- Alarm setting for all parameters

Select alarm key -> Alarm Settings

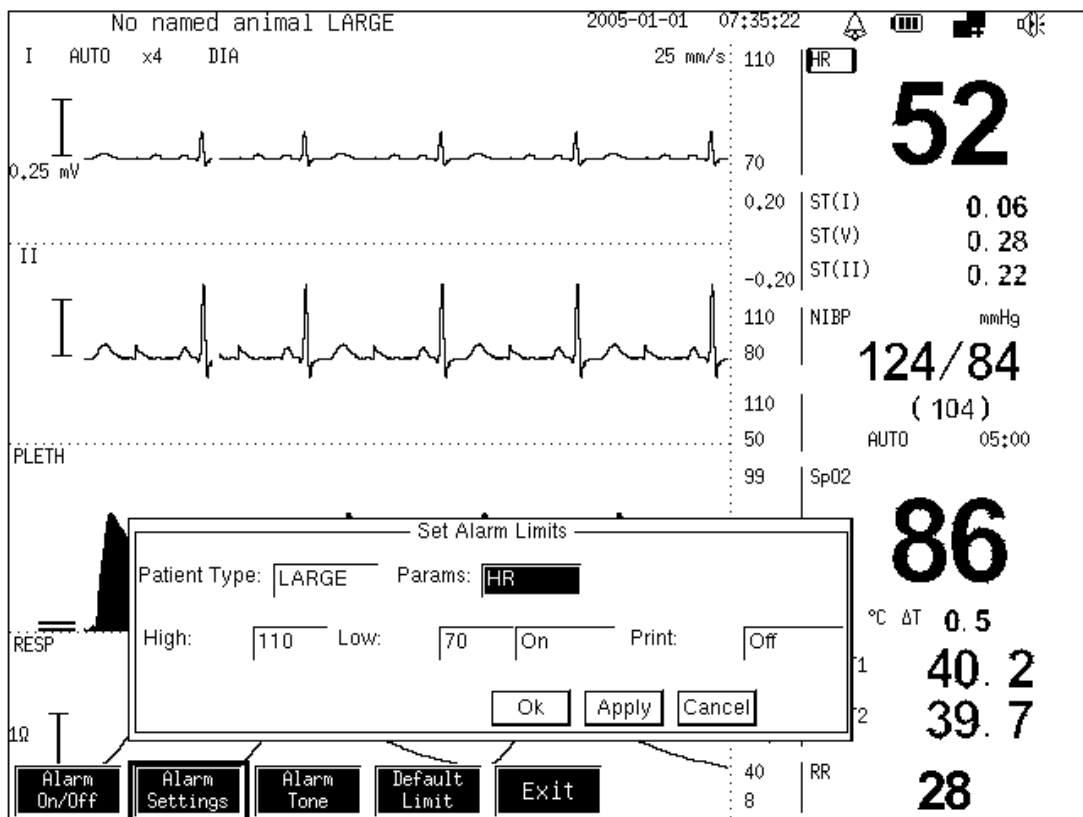

3.10 Preparatory work

- Check the Monitor

Before the measurement, please check whether the Monitor has mechanical failure, which including all external cables, plug-ins and accessories connecting to it.

- Connecting the power

After connecting the power, press the “on/off” Key on the Monitor.

- Setting modules

Configure the modules which ought to be observed, check whether the animal cables and sensors are properly connected.

- Starting monitoring


Connect the animal to the Monitor. Verify the “Monitor setting”, “Alarm limit”, alarm volume, animal type and pacing impulse status to be suitable to you animal, if necessary, change them. For the detailed information for measuring your required parameters, please refer to the chapters and sections about the measurement.

- Cut off the power

Press the “on/off” Key, and then pull off the power cord.

3.11 Remote maintenance

Select menu key -> System Setup -> Maintenance -> Remote maintenance

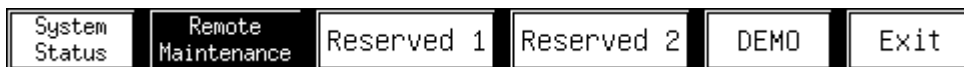


Enter the Remote Maintenance status, the Central Management System can remotely upgrade the Monitor.

3.11.1 Setting your monitor

Menu settings

■ Menu

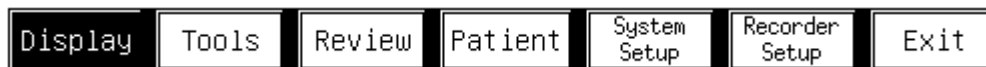

<Display> mainly set the display interface.

<Tools> mainly set the practical tools.

<Review> mainly observe the historic data.

<Animal> mainly set animal's information.

<System Setup> mainly set the system information.


<Recorder Setup> mainly set the printer parameters.

■ Display menu




<Fixed Format> it can set the display channel as the 1 ECG, 2ECGs, 3 ECGs, 7ECGs, Huge Digit.


- 1 ECG

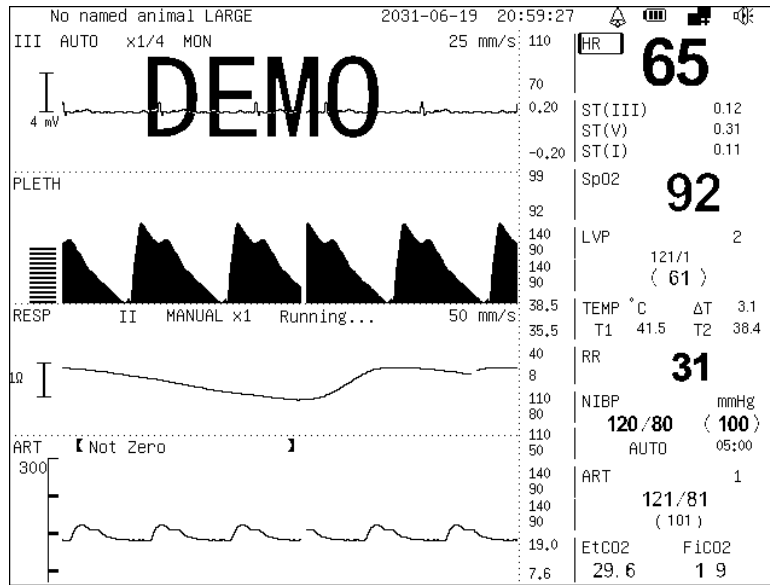


- 2 ECGs

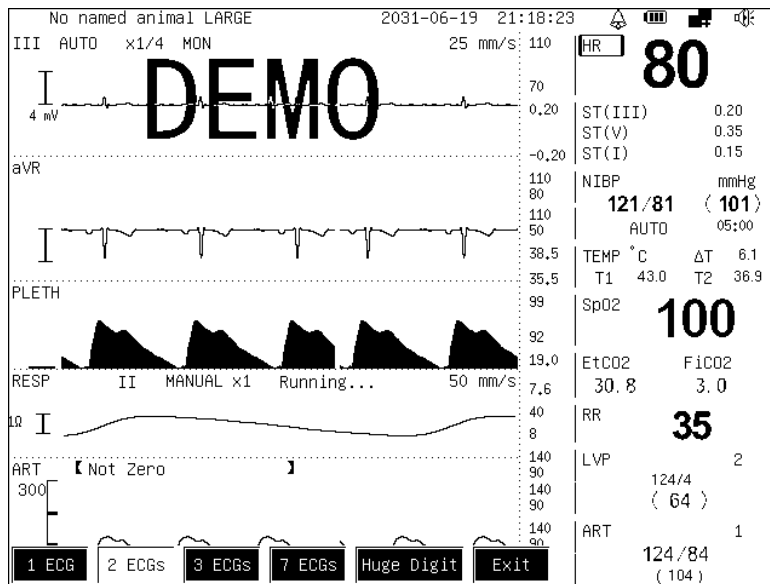


- 3 ECGs

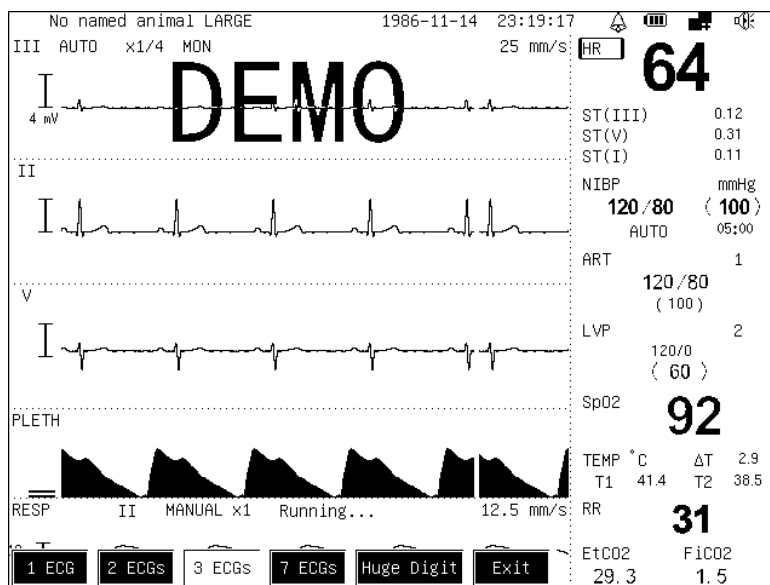

- 7 ECGs

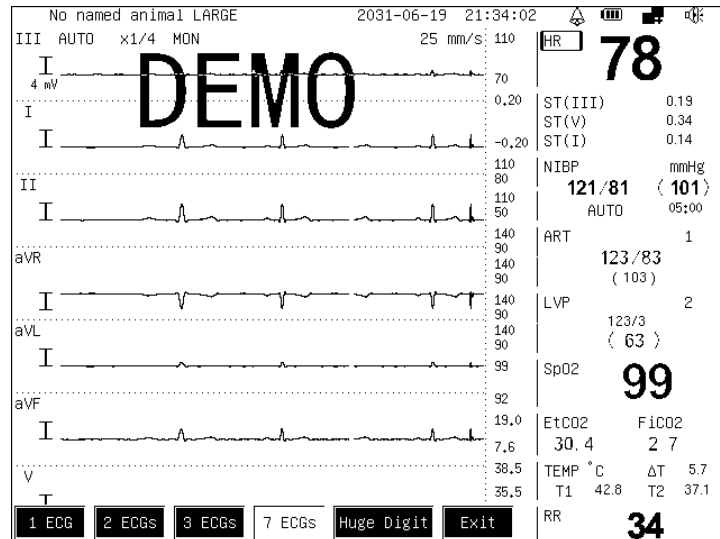


- Huge Digit

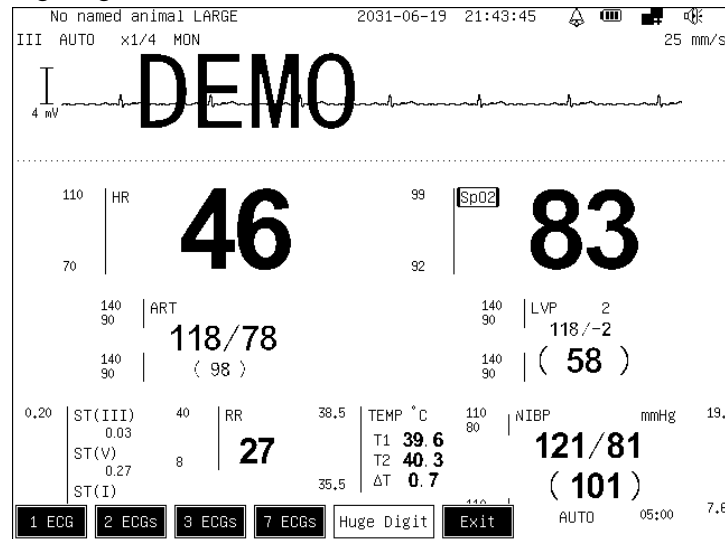

<User Format> it can set the display channel as the channel format defined by the user.


It supports five user-defined channel formats.

<Format Setup> the user can define the display channel by himself.

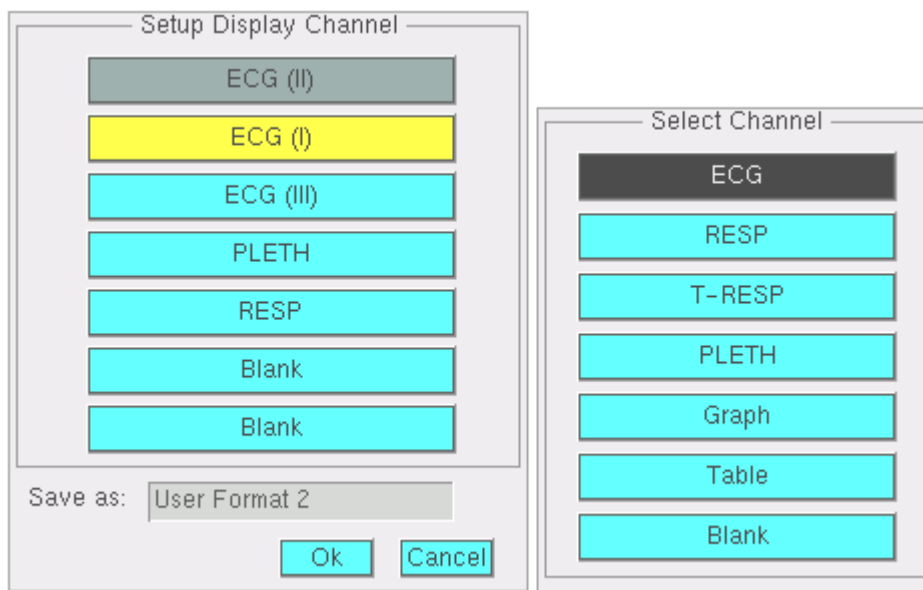


- Note:
- The same channel cannot be displayed repeatedly.
 - Blank means no channel to be displayed at this position.
 - Select the ECG channel, the system will automatically configure the lead of the ECG channel.
 - The configured channel can be saved as the user-defined format, and later this format can be selected at any time from the <User Format>.


<oxyCRG> display the oxyCRG chart.

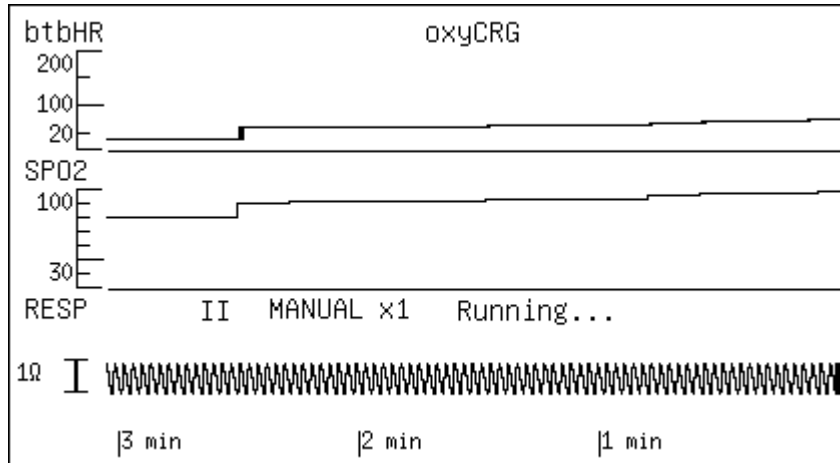

<Alarm Limit> display or hide the alarm limit parameter area.

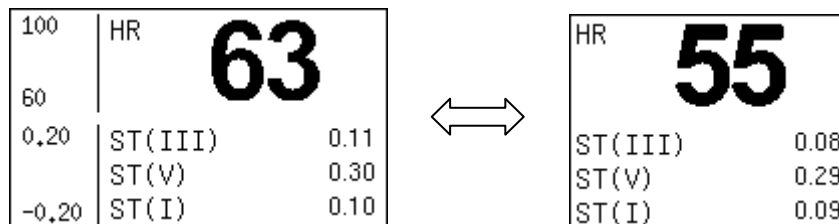


<Waveform Speed> it can adjust the waveform speed of the ECG / PLETH/IBP and RESP/CO2.

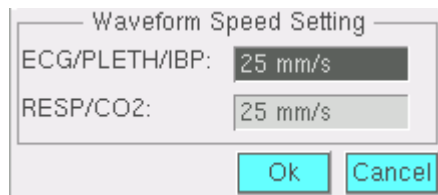

<Other Settings> it can set other display formats.

- NIBP Size
Three display formats are switchable.

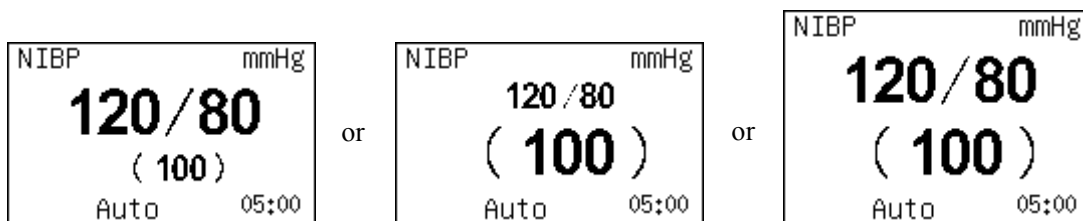

- IBP1 Size

Three display formats are switchable.




- IBP2 Size

Three display formats are switchable.




- IBP Overlap

Display the two IBP waves in a coordinary.

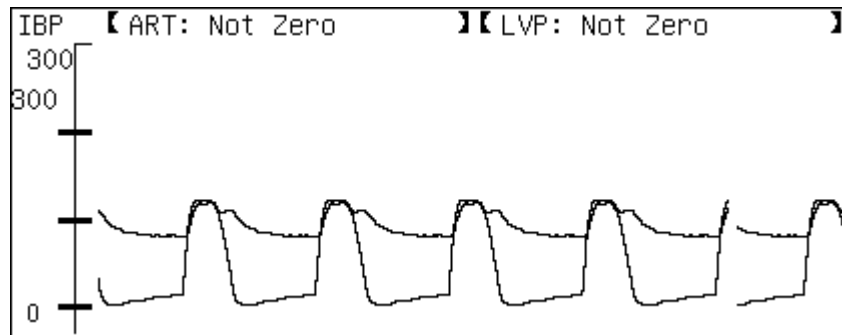

- Menu Font

Menu font can be switched to large or small font.

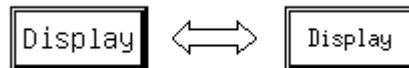


- Screen Font

The screen font can be switched large or small font.




- Pelth Fill

Pelth waveform can be switched to filling or blank.

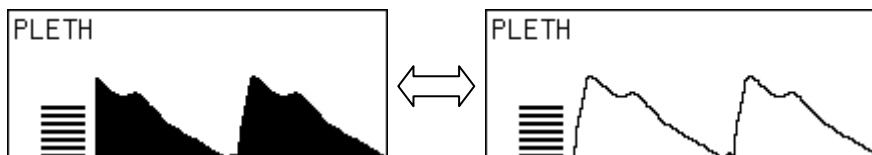

■ Tools menu

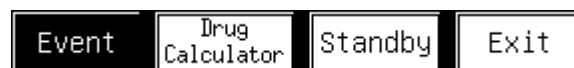

<Event> it can set the events.



<Drug Calculator> the drug calculator can be called.

<Standby> it makes the Monitor into standby mode.

■ Review menu

<Trend> you can observe the trend data, for more detailed information, refer to Chapter of Trend Analysis.

<Recall> it can replay the waveform, for more detailed information, refer to Chapter of the Waveform Freeze and Replay.

■ Animal menu


The animal information can be set.

■ System Setup menu

<Set Time> you can set the current time for the Monitor.


<Region Settings> you can set the parameters relating to the regions.

- Select Units

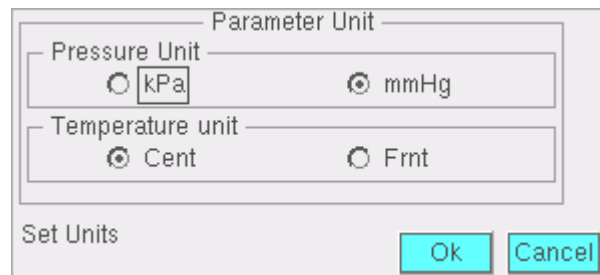


The dialog box is titled "Parameter Unit". It contains two sections: "Pressure Unit" with radio buttons for "kPa" and "mmHg" (selected), and "Temperature unit" with radio buttons for "Cent" (selected) and "Frnt". At the bottom, there are "Set Units", "Ok", and "Cancel" buttons.

- Date Format

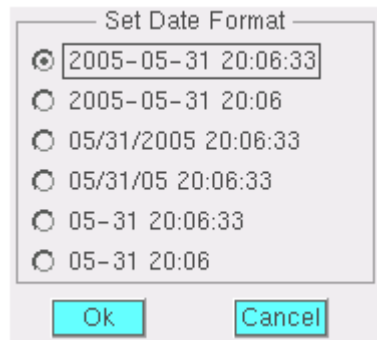

The dialog box is titled "Set Date Format". It lists seven date and time format options with radio buttons. The first option, "2005-05-31 20:06:33", is selected. At the bottom, there are "Ok" and "Cancel" buttons.

- Power Frequency

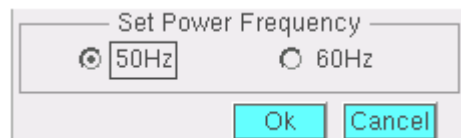


The dialog box is titled "Set Power Frequency". It has two radio button options: "50Hz" (selected) and "60Hz". At the bottom, there are "Ok" and "Cancel" buttons.

Note: please set the right power frequency, or the signal noise of the ECG will be loud.

- Language

It can be set in Chinese or English.


<Network Settings>

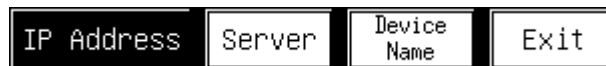

A horizontal menu with four buttons: "IP Address", "Server", "Device Name", and "Exit".

- IP Address

You can set the IP address of the Monitor.

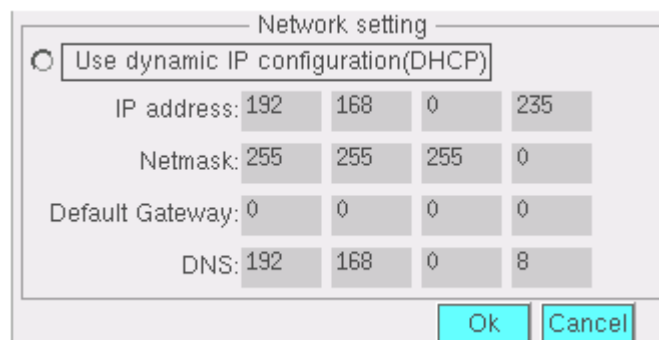

The dialog box is titled "Network setting". It has a radio button for "Use dynamic IP configuration(DHCP)" which is selected. Below this are four rows of input fields for IP address, Netmask, Default Gateway, and DNS, each with four numerical boxes. The values are: IP address: 192, 168, 0, 235; Netmask: 255, 255, 255, 0; Default Gateway: 0, 0, 0, 0; DNS: 192, 168, 0, 8. At the bottom, there are "Ok" and "Cancel" buttons.

- Server

You can set the IP address and domain name of the Server.

- Device Name

You can set the name of the Monitor.

<Calibration> you can set the calibration function.

<Maintenance> you can set the maintenance function.

Caution: This maintenance function only can be set by the Administrator.

The system maintenance function includes: System Status, Remote Maintenance, Reserved1, Reserved 2 and DEMO.

The maintenance menu is as following:

The System Status: to check the relevant system status. The menu is as following:

[Main Settings] list sample is as following:

[Alarm Limits] list sample is as following:

Main Limits				ST Alarm Limits		
Param	On/Off	Low limit	High limit	Lead type	Low limit	High limit
HR	Off	70	110	I	-0.22	0.22
SPO2	On	92	99	II	-0.20	0.20
RESP	On	8	40	III	-0.20	0.20
TEMP1	On	35.5	38.5	aVR	-0.22	0.22
TEMP2	On	35.5	38.5	aVL	-0.20	0.20
SBP	On	80	110	aVF	-0.20	0.20
MAP	On	70	125	V	-0.20	0.20
DBP	On	50	110			

Ok

[Software Version] list sample is as following:

Software Information
Serial No: G3F5050075
GBA440: GMC13D01.68
GAA718: GM033000.06
KERNEL: 2.4.18-1.01
Mon, 14 May 2007 15:35:40 +0800
Digital Oximeter: yes
Digital Oximeter Software Version: 1.5
COMM: 262690,130,0,9878
CO2: 0,0,0,0
No:

Ok


You shall pay more attention on the following operation functions:

Remote maintenance: it is used to upgrade the Monitor software;


Reserved 1, Reserved 2;

Demo Switch: to turn on the demonstration status.

<Default Settings> the ex-works settings can be called.

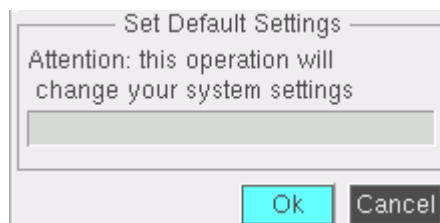

■ Alarm Setup menu

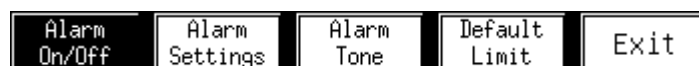

<Alarm On/Off> open or close the alarm of all the parameters, the detailed information shows in 3.9 Setting Alarm.


<Alarm Settings> set the alarm limits of the parameters, the detailed information shows in 3.9 Setting Alarm.

<Alarm Tone> set the pulse and alarm tone, shows as following:



The screenshot shows a dialog box titled "Set Pulse & Alarm Tone". It contains two rows of controls. The first row is for "Pulse Tone", with a numeric input field containing "800" and a cyan "Test" button to its right. The second row is for "Alarm Tone", with a numeric input field containing "900" and a cyan "Test" button to its right. At the bottom right of the dialog, there are two cyan buttons: "Ok" and "Cancel".

<Default Limit> can set the value of the alarm parameters as the default value, shows as following:

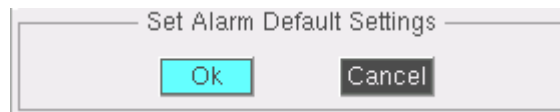

The screenshot shows a dialog box titled "Set Alarm Default Settings". It contains two cyan buttons: "Ok" and "Cancel".

Chapter 4 ECG Monitoring

This Monitor adopts 5-lead ECG cable to collect the ECG signal from the animal, and it can simultaneously display seven leads of ECG signal, ECG channel displayed on the top part of the screen is the major channel of the ECG signal, the lead type selected by this Monitor we called it the major lead, the Monitor calculates the heart rate value according to the data of this lead, it also controls the gain (if set in automatic mode) according to the lead data. All leads adopt the same gain control and same measurement mode, the Monitor can separately check whether the connections of the leads fall off, and display them on the screen.

4.1 ECG electrode monitoring

In order to monitor the animal's ECG, we use a five-electrode ECG cable by which we can get more than 12 selective ECG lead, while the certain lead of the Monitor is effective, the animal's ECG waveform corresponding to the lead will be displayed.

The ECG cable includes two parts: the main cable connecting the Monitor and the electrode lead line connecting the animal. Usually the ECG cable has five electrodes. Pay attention to the following items when installing the electrodes:

1) The selection and use of the electrode

Generally the electrode for monitoring is disposable electrode made from Ag-AgCl (silver - silver chloride). But before use, you shall confirm that the electrode is within the valid date. If the unqualified electrode is used, it will lead to inaccurate monitoring result.

2) Skin pretreatment

Dogs and cats have an adequate amount of electrolyte material in their skin and hair, so moistening lead sites with isopropyl alcohol will be enough, although it is suggested that before installing the electrode, you shall carry out skin pretreatment for animal at first.

- (1) Wet the hair and skin at the lead attachment site with isopropyl alcohol.
- (2) Install the electrode on the animal. If monitoring for more than 30' to 60' the use of electrical ointment is advisable.
- (3) For animals with dense undercoat, rub ointment with fingers to assure that it has made contact with skin.
- (4) Connect the electrode lead with the animal cable.
- (5) Crocodile clips are included with the equipment. Be sure to open them wide enough to grasp the skin. Take care not to cause injuries.
- (6) Confirm the Monitor is power on.

4.1.1 Installing the ECG electrode

For the installation of the five-lead ECG electrode, refer to the ECG electrode installation figure as shown in Figure 4.1.

Note: the sign name and color of the electrode is American standard, while inside the square brackets (-) are European standard.

The electrode installation for five-lead device

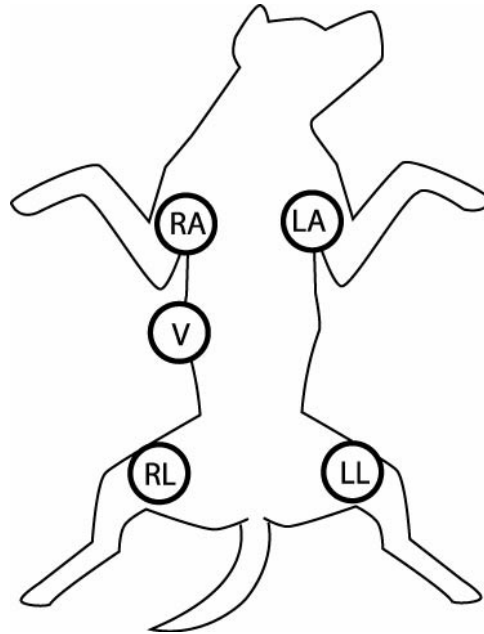


Figure 4.1 Five-lead ECG electrode installation positions

- White (red) electrode (RA) — Right foreleg
- Black (yellow) electrode (LA) — Left foreleg
- Green (black) electrode (RL) — Right hind leg
- Red (green) electrode (LL) — Left hind leg
- Brown electrode (V) — Fourth intercostal space


Notice

- Select the place with stable ECG signal or the little inference of skeleton activities to place the electrodes.
- The non-grounding equipment near the animal and the electrosurgery unit (ESU) interference may lead to waveform problems.
- When using the surgical electrocauterize unit, it prohibited to put the electrode near the grounding of the surgical electrocauterize unit, or a lot of interference on the ECG signal will emerge, the electrode spaces shall be equal in axes to the electrocauterize blade and electrocauterize grounding.
- The electrode shall be fixed properly to ensure the reliable contact of the conductor with skin.
- In order to ensure the animal's safety, all leads shall be connected to the animal.

4.1.2 ECG monitoring interface

4.1.2.1 Introduction on the main lead

On the Monitor interface, the ECG waveform is green, including two kinds of ECG channel, of which, the top ECG channel is the main channel of the ECG, we called it the main lead channel, as shown in following Figure:

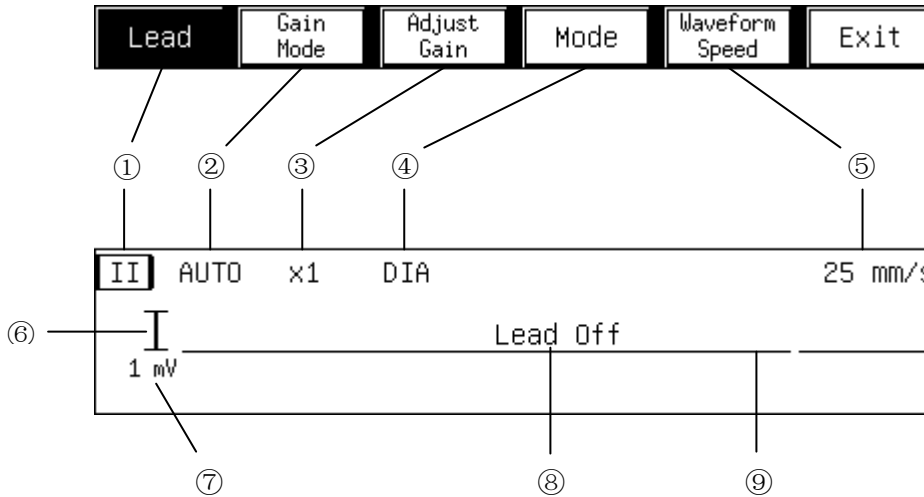

- ①main lead type ②gain mode ③gain ④measurement mode
 ⑤waveform speed ⑥scale ⑦scale range ⑧lead off sign ⑨ECG waveform

Select this ① position with the cursor, press the rotating mouse and the ECG menu shows. The function of the knobs on the ECG menu are introduced here:

【Lead】 : Switch the lead type. Press the mouse once, it switches once. There are seven lead types such as I, II, III, aVR, aVL, aVF, V for selection. The main lead is a very important lead, in order to get accurate heart rate and gain control performance, the doctor shall select the lead of large range and little noise as the main lead.

Notice: in order to avoid emerging leads repeatedly, when switching the lead type in the current channel, the lead types in the other channels will be changed automatically.

【Gain Mode】 : Switch to automatic gain and manual gain, press the mouse once, it switches once.

Gain means the magnification for the ECG waveform signal, this Monitor have five magnification, they are x1/4, x1/2, 1, 2, 4. the x1 in the center is one time magnification, under this magnification, the range ⑦ of the scale⑥ at the left side of the ECG waveform is 1mV. Under the standard graph display format, the length of the scale⑥ is duly 10mm, then the waveform of the 10mm displayed on the standard interface means the ECG signal of 1mV. Then under the condition of x2 magnification, the waveform of the 10mm displayed on the standard interface means the ECG signal of 0.5 mV, and so on, for detailed data refer to the following table:

Gain factor	Automatic gain	Manual gain	10mmscale range
X1/4	Adopted	Applicable	4 mV
X1/2	Adopted	Applicable	2 mV
X1	Adopted	Applicable	1 mV
X2	Adopted	Applicable	0.5 mV
X4	Adopted	Applicable	0.25 mV

Notice: All leads including the main lead and other leads adopt the same gain control.

This Monitor provides two methods to regulate the electrocardiogram waveform range:

Automatic mode: Under the automatic mode, according to the main lead waveform data, the Monitor will automatically regulate the gain to amplify the ECG waveform of the main lead without the distortion as possible. The defect of the automatic mode is slow regulation speed, under this mode, the manual gain regulation knob **【Adjust Gain】** is voidable.

Manual mode: Under the manual mode, the Monitor will not automatically regulate the ECG gain value. At this time, the ECG gain is regulated via the **【Adjust Gain】** knob. The feature of the mode is fast response, the waveform will change immediately after the gain is regulated, and can amplify the waveform to a distortion status to display.

These two modes are switched via the knob **【Gain Mode】**.

【Adjust Gain】 : Manually regulate the ECG gain. After pressing down the rotating mouse, regulate the gain via the rotating mouse. Rotate to left, the gain and wave range become small, when rotate to right, the gain and wave range become large, after accomplishment of the regulation, press down the rotating mouse again.


【Mode】 : The measurement mode means the filter bandwidth of the ECG channel. Under the diagnosis mode, the filter bandwidth is set to the widest, which can get the detailed information of the ECG signal. The more detailed information is helpful to more accurately judge the status of the ECG signal. But when the importation of the detailed information, the environmental noise may be imported tot, one example for this kind of noise the HF electrotome, the noise mixes with the real ECG signal, so that the ECG signal can not be distinguished. In order to adapt to the noise interference situation, the Monitor also provides two measurement modes, namely the monitoring mode and the operation mode, under these two modes, it adopts narrower bandwidth to measure and gains more smooth signal, the user can make choices according to the exact situations, the comparing data are as followings:

Measurement mode	Bandwidth	Details	Noise
Diagnosis	0.05Hz~100Hz	Most	Most
Monitoring	0.5Hz~40Hz	Medium	Medium
Operation	1.0Hz~25Hz	Least	Least

Notice: all leads including the main lead and other lead adopt the same measurement mode.

【Waveform Speed】 : regulate the waveform speed for parameters relating to the heart, including all ECG channels and blood volume waveform channels, it has four speeds such as 6.25 mm/s, 12.5 mm/s, 25 mm/s, 50 mm/s.

4.1.2.2 Introduction on the interface of other lead channels

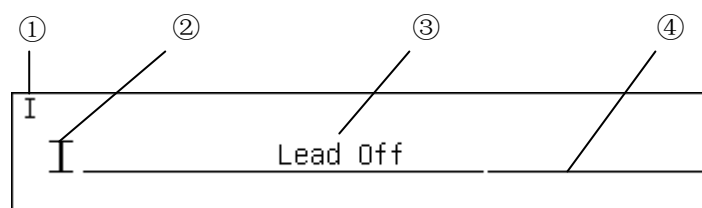


①lead type ②scale ③lead off sign ④ECG waveform

Select the position ① with the cursor, press the rotating mouse and switch lead types.

4.2 Parameter display

ECG parameters include the heart rate and ST segment, except the large character mode, the heart rate and ST segment are displayed on the right corner of the screen, the format is as shown in the Figure below.

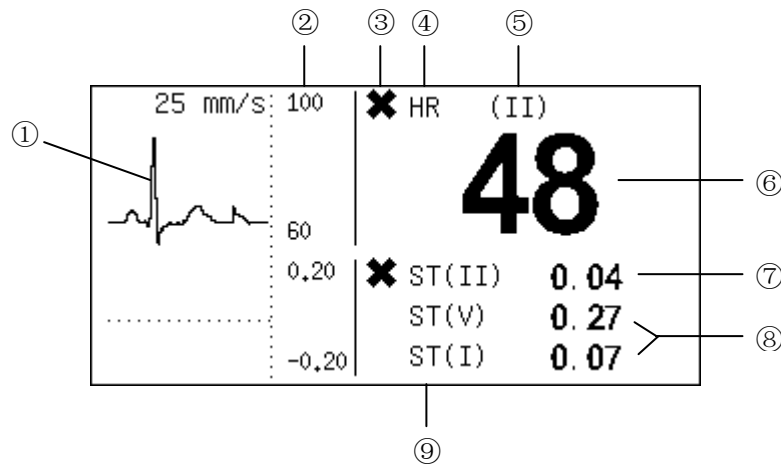


①ECG waveform ②heart rate and main lead ST alarm limit ③heart rate alarm off sign ④heart rate label ⑤heart rate source sign ⑥heart rate ⑦main lead ST segment ⑧other lead ST segment ⑨ST segment label. Tune the rotating mouse, move the cursor, select the heart rate label, and press down the rotating mouse, pop up the heart rate setting menu, as shown in the Figure below:




【Alarm On/Off:】 : Switch between the alarm on and alarm off.

【Alarm Setting】 : Check and regulate the alarm setting of the heart rate.


【HR Source】 : Set up the selection strategy for the heart rate source. Press down and the following menu shows:

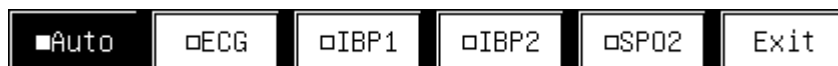

【AVRG.】 : Set up the calculation average periods of the heart rate and ST segment, Press down and the following menu shows:

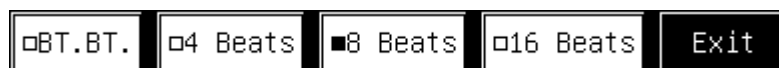

【ST On/Off】 : Open and close the ST segment display, as shown in the Figure below:

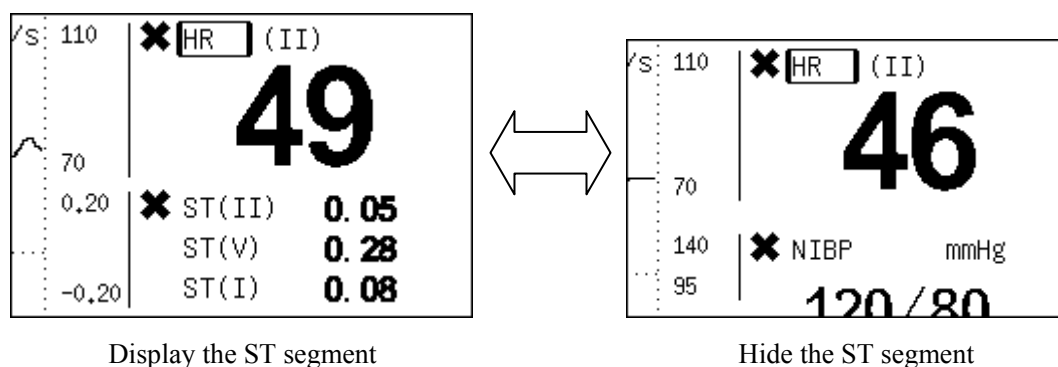

Heart rate source selection function:

During the oxygen saturation testing process, the Monitor can measure the pulse rate. Pulse rate and heart rate, one of them can be displayed on the position of the heart rate, and this selection can be called heart rate source selection. This Monitor supports two kinds of heart rate source selection system such as automatic selection and manual selection, the current selected value is displayed here, the heart rate source mark and heart rate label clearly show the current heart rate source, the list is as following:

Label	The color of the heart rate value	Heart rate source mark	Source
Heart rate	Green	I, II, III, aVR, aVL, aVF, V	Form ECG signal
Pulse rate	Red	SpO2	From SpO2 signal

Notice: pulse rate and heart rate use one set of alarm limit and alarm on/off. In spite where the heart rate value comes from, the alarm systems adopted are totally the same.

4.3 Maintenance

Caution If the ECG cable is damaged or aging, it shall be replaced by a new one.

4.4 Trouble shooting

Problem	Possible reasons	Suggestions and way-out
Lead off	The ECG electrode fallen off from the animal's body or the ECG cable fallen off from the Monitor.	Ensure the electrode, lead and cable are properly connected.
ECG signal noise is too loud.	The measurement mode is not correctly set according to the environment. The power frequency is not set in accordance with the local power frequency. The Monitor is bad grounding. The electrode is badly connected to the animal.	Change the electrocardiogram measurement mode. Set the power frequency right, refer to Chapter Three. Check the grounding system of the Monitor. Keep the animal quiet, and guarantee the reliable connection of the electrode.
No heart rate display No heart beat sound	The ECG signal is very weak. (<0.25mV) The electrode is badly connected to the animal. The animal type is not correctly set. The gain is not enough.	Check the connection between the electrode and animal, ensure the reliable touch of the electrode and animal's skin. Correctly set the animal type. If it is the manual gain regulation, please tune the gain to the right magnification.

Chapter 5 NIBP Monitoring

5.1 Overview

- The noninvasive blood pressure (NIBP) measurement adopts the auscultation.
- It is applicable to dogs and cats and other animals.
- The measurement mode: manual measurement, automatic measurement and consecutive measurement, each mode displays the systolic blood pressure, mean pressure and diastolic blood pressure.

5.2 Blood pressure cuff installation

Make sure that the blood pressure cuff has been completely deflated; fix the blood pressure cuff on the upper arm or thigh of the animal as shown in Figure 5.1.

Ensure the pipe exit of the blood pressure cuff is rightly on the proper artery, confirm the degree of tightness of blood pressure cuff, don't be too tight, or it can result in discoloration or ischemia of the extremity.

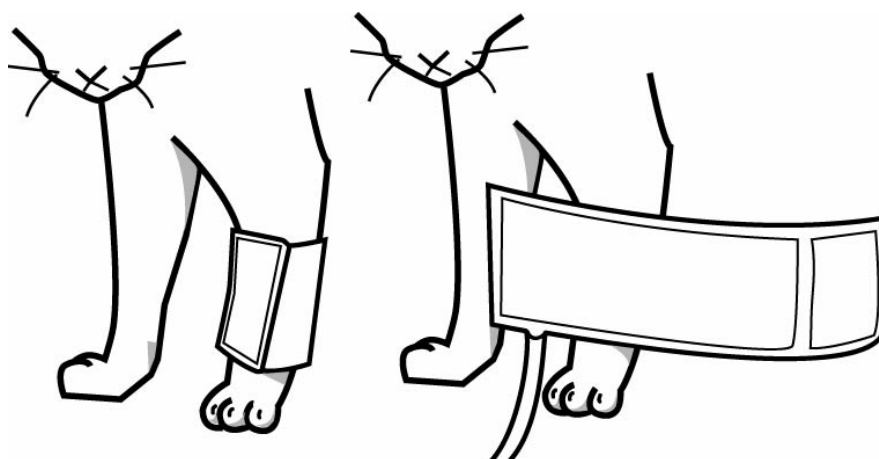


Figure 5.1 blood pressure cuff installation position

Cats must be left in its owners lap to keep them calm. Try to have the animal in an area away from noise and bright lights. Front limbs are preferable sites although in some nervous animals the tail could be an alternative location.

Sternal, lateral or dorsal recumbent positions are the choice for NIBP measurements in dogs although this could be a position difficult to achieve in large dogs which don't cooperate. If so, place the dog in sitting position and take measurements from the metacarpus of any of the front limbs.

Other areas in dogs could be the metacarpus and the anterior tibia or the tail. Clip hair if it is too thick and could affect measurements. Avoid using the tail during anesthesia.

Notice:

- ① The width of the blood pressure cuff should span 40% of limb circumference or two-thirds of the upper arm length, the section for air charging shall be long enough to round 50~80% of the arm or leg, the inappropriate size of the blood pressure cuff will lead to error reading.

- ② The limb to be measured shall be in the same level with the heart, if it can not be realized, the following correction shall be carried out on the measurement result.
- If the blood pressure cuff is above the heart level, add 0.9mmHg (0.10kPa) to the displayed value per centimeter, or 2.2 mmHg (0.25kPa) per inch.
 - If the blood pressure cuff is low to the heart level, subtract 0.9mmHg (0.10kPa) from the displayed value per centimeter, or 2.2mmHg (0.25kPa) per inch
 - If you have suspicious on the accuracy of the reading, before checking the function of the Monitor, check the animal's vital signs with possible means at first.

Caution

- ⚠ Please use the accessories approved by our company, or it will lead to the malfunction of the Monitor.
- ⚠ Do not use it to measure the NIBP on the animal with skin damage.
- ⚠ For those animals with severe blood coagulation mechanism malfunction, the automatic blood pressure measurement shall be adopted or not on the basis of the clinic evaluation, because the limb friction with the blood pressure cuff may lead to the risk of haematoma.
- ⚠ Do not install the blood pressure cuff on the limb with intravenous transfusion or pipe, because when the blood pressure cuff is charged with air, it will lead to the organ damage around the pipe due to slow intravenous transfusion or block.
- ⚠ Under the automatic mode, if the NIBP measurement time is too long, the friction between the blood pressure cuff and limb will lead to the purpura, lack of blood or nerve damage. When monitoring the animal, the veterinarian shall often check if the color of the remote limb end is normal, the warm degree and sensitivity. If some abnormal situations occur, stop the blood measurement.
- ⚠ The air charging pipe between the blood pressure cuff and the Monitor shall be unobstructed and cannot be tangled.

5.3 The limitations on the measurement

The blood pressure measurement of this Monitor is vibration measurement. This measurement looks for a regular artery pressure pulsation. If animal's condition is not applicable to use this method, the measurement value may become unreliable; the measurement time will be increased too. The user shall realize that the following situation interfere the measurement method, which make the measurement becomes unreliable or the longer measurement time, under this condition, the animal's status will make the measurement impossible.

Animal's exercises: for example, if the animal is exercise, shivering or jerking, which will make the measurement result unreliable or even impossible, these conditions will interfere in the artery pressure pulsation measurement, further more, the measurement time will be prolonged.

Arrhythmia: for example, if the animal is in irregular heart beating due to the cardiac arrhythmia, the measurement result will be unreliable or the measurement can not be carried out, the measurement time will be prolonged too.

Pressure variations: if the animal's blood pressure changes rapidly during certain measurement process, the measurement will be unreliable or can not be conducted,

Severe shock: if the animal is in severe shock or very low body temperature, the measurement will be unreliable, because the blood flowed to the periphery, which leads to the decrease of the artery pulsation.

Ultimate heart rate: when the heart rate is low than 30bpm (heart beat /minute) or is more than 300bpm, the measurement cannot be carried out.

5.4 Normal Blood pressure values

Feline normal values

The blood pressure values for cats are the same for all breeds.


Normal feline blood pressure: 124/84.

Canine normal values

The average canine blood pressure is: 133/75, although these seem to be breeds specific. The following table shows normal values for a certain number of breeds.


Breed	Systolic (mmHg)	Diastolic (mmHg)	Pulse Rate
Labrador Retriever	118 ± 17	66 ± 13	99 ± 19
Golden Retriever	122 ± 14	70 ± 11	95 ± 15
Great Pyrenees	120 ± 16	66 ± 6	95 ± 15
Yorkshire Terrier	121 ± 12	69 ± 13	120 ± 14
West Highland	126 ± 6	83 ± 7	112 ± 13
Border Collie	131 ± 14	75 ± 12	101 ± 21
King Charles Spaniel	131 ± 16	72 ± 14	124 ± 24
German Shepherd	132 ± 13	75 ± 10	108 ± 23
Terrier	136 ± 16	76 ± 12	104 ± 16
Bullterrier	134 ± 12	77 ± 17	122 ± 6
Chihuahua	134 ± 9	84 ± 12	109 ± 12
Miniature Breeds	136 ± 13	74 ± 17	117 ± 13
Pomeranian	136 ± 12	76 ± 13	131 ± 14
Beagle	140 ± 15	79 ± 13	104 ± 16
Dachshound	142 ± 10	85 ± 15	98 ± 17
Saluki	143 ± 16	88 ± 10	98 ± 22
Greyhound	149 ± 20	87 ± 16	114 ± 28
Pointer	145 ± 17	83 ± 15	102 ± 14

5.5 Monitoring display interface

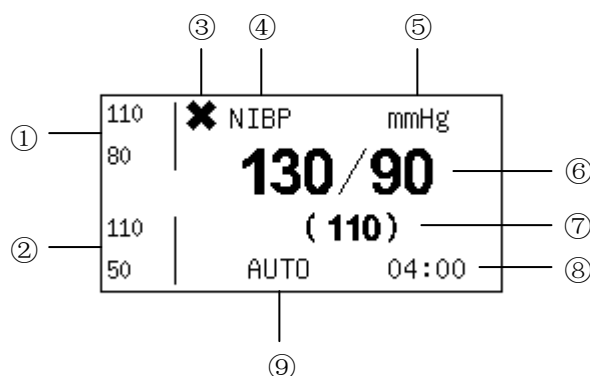

- ① Systolic blood pressure alarm limit.
- ② Diastolic blood pressure alarm limit.
- ③ NIBP alarm off sign: when the alarm for systolic blood pressure, diastolic blood pressure and mean pressure is in prohibition status, this sign will be displayed, if one parameter is allowed to alarm, then this sign will disappear.
- ④ NIBP, press down and NIBP setting menu shows.
- ⑤ Pressure unit.
- ⑥ Systolic blood pressure and diastolic blood pressure.
- ⑦ Mean pressure.
- ⑧ ⑨ The information prompting bar for the NIBP. Different information is displayed on it according to different situations.

The NIBP setting menu is as shown in the Figure below:

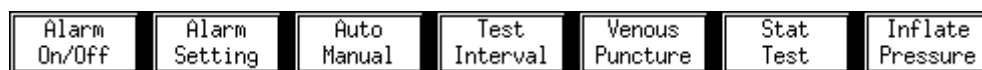

- 【Alarm On/Off】 simultaneously switch the alarm on/off for systolic blood pressure, diastolic blood pressure and mean pressure.
- 【Alarm Setting】 check and set the alarm for systolic blood pressure, diastolic blood pressure and mean pressure.
- 【Auto/Manual】 switch the automatic measurement mode and manual measurement mode for the noninvasive blood pressure
- 【Test Interval】 set the automatic measurement time span under the automatic measurement mode
- 【Venous Puncture】 enter and exit the venipuncture status.
- 【Stat Test】 start up the continuous measurement.
- 【Inflate Pressure】 set the initial air charging pressure for the manual measurement mode.

5.6 Introduction on functions

The noninvasive blood pressure has two functions, one is blood pressure measurement, and the other is the venipuncture function.

5.6.1 Blood pressure measurement

When the blood pressure measurement process starts up, the Monitor charges the air into the blood pressure cuff, the human body blood pressure, including the systolic blood pressure, diastolic blood pressure and mean

pressure are measured out via the sensors, a measurement process is about 40 seconds.

There are three start-up methods for blood pressure measurement, one is manual measurement mode, the other is automatic measurement mode, and the third is consecutive measurement mode.

- **Manual measurement mode:** under the manual measurement mode, each blood pressure measurement process shall be started up manually; the start up method is pressing the <START> on the panel of the Monitor, under this mode, the blood pressure cuff air charging pressure is set manually.

The manual measurement mode is entered by using **【Auto/Manual】** knob. When the "MANUAL" displays in the position ⑨, it means that it enters into the manual measurement mode. Position ⑧ does not display any information.

- **Automatic measurement mode:** under the automatic measurement mode, the Monitor will periodically start up the blood pressure measurement according to the assigned time span. If you want to enter the automatic measurement mode, use the **【Auto/Manual】** knob to switch, when position ⑨ displays "AUTO", it means it has entered the automatic measurement mode, now position ⑧ shows clock format of the automatic measurement time span.

The automatic measurement mode has two status, one is the pause status, at this time, the clock at position ⑧ stopped, it means that although it has entered the automatic mode, but the automatic measurement cycle has not been started yet; the other status is the operational status, at this time the clock at position ⑧ runs, it means the automatic measurement cycle has been started. The clock at position ⑧ indicates how much time still remains for next measurement. In the pause status, carry out one manual blood pressure measurement, then it enters the automatic measurement operational status, you can see that clock at position ⑧ begins to run, if you want to stop the automatic measurement status, you can switch to the manual mode, then switch back to the automatic mode. If you press down the <START> knob during the blood pressure measurement, it not only can stop the measurement immediately, but also can switch the automatic measurement operational status to the pause status.

Notice: whether you are in the automatic mode or manual mode, press the <START> knob on the panel, you can start up the blood pressure measurement.

Notice: during the blood pressure measurement process, press down <START> knob, the Monitor stops the blood pressure measurement immediately and deflate the blood pressure cuff at the same time.


- **Continuous blood pressure measurement**

The blood pressure continuous fast measurement mode is a mode used for operation or emergent animal, its feature is that if it is started, it will carry out blood pressure measurements time after time for consecutive five minutes, under this mode, the measurement mode adopted by the Monitor is simpler than other modes, after five minutes, the Monitor restores to the status (automatic or manual) before the continuous fast measurement.

5.6.2 Venipuncture

Under the venipuncture mode, the Monitor charge air into the blood pressure cuff to the assigned pressure, it keeps the assigned pressure during the assigned time length, then deflate the air.

The setting interface for the venipuncture is as following:

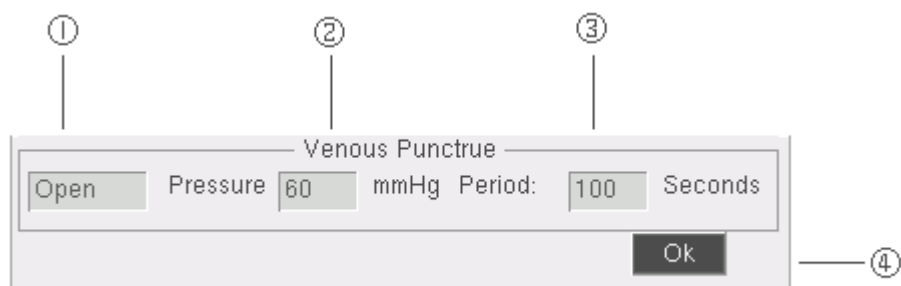


- ① The venipuncture mode switch.
- ② Keeping pressure.
- ③ Keeping time. The keeping time is calculated from the air charging time, until the air is deflated.

Animal type	Pressure setting range	Default pressure	Maximum continuous time
Large	20-80 mmHg	40 mmHg	100 s
Small	20-50 mmHg	30 mmHg	60 s

- ④ After setting, press this key.

When entering the venipuncture mode, the prompting bar will prompt the information as shown in the Figure below, it means it enters the venipuncture mode:

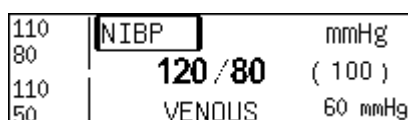

Under the venipuncture mode, you can use the <START> knob to start and stop the air charging.

Notice: the blood pressure measurement has close relationship with the animal type, please set correct animal type to ensure safety, refer to the relevant content concerning the Monitor setting.

5.7 Maintenance and cleaning

■ The reusable blood pressure cuff

The blood pressure cuff can be put in the ordinary hot air oven with high pressure vapor, it can be sterilized with gas or radiation, or put into the decontamination liquid to kill the bacteria, but remember, when you use this method, take out the rubber bag please. The blood pressure cuff cannot be dry-cleaned. The blood pressure cuff can be washed with hand or machine, its life time can be prolonged by hand washing. When you washed it, after the blood pressure cuff is completely dry, put the rubber bad into it.

If you want to put the rubber bag into the blood pressure cuff again, at first you can put the rubber bag at the front end of the blood pressure cuff, thus, the rubber pipe and the large slot of the long end of the blood pressure cuff are in the same line, now, roll the rubber bag longitudinally and insert into the large slot of the blood pressure cuff, hold the rubber pipe and blood pressure cuff, shake the whole blood pressure cuff until the rubber bag get into its position, then take the rubber pipe into the blood pressure cuff, and get it out via the inner of the hole.

■ Disposable blood pressure cuff

The disposable blood pressure cuff is appointed to be used on one animal, do not use the same blood pressure cuff on other animal, the disposable blood pressure cuff can not be sterilized or kill the bacteria with high

pressure vapor, the disposable blood pressure cuff can be washed with soap and clean water to control the infection.

Caution

- ⚠ Do not press or limit the rubber pipe of the blood pressure cuff.
- ⚠ Do not allow the water or detergent liquid get inside the air passage connection on the front of the panel, or the instrument may be damaged.
- ⚠ When the reusable blood pressure cuff is not connected to the instrument or being washed, be cautious to prevent the liquid from entering the rubber pipe and being soaked into the instrument.
- ⚠ When clean the Monitor, only clean the external of the connector socket, and do not clean its internal part.

Note:

- ⚠ In order to protect the environment, the disposable blood pressure cuff shall be recycled or be disposed properly.

5.8 Trouble shooting

Problems	Possible reasons	Suggestions and way-out
Start failure	The Monitor hardware error	Stop the NIBP measurement function, inform the supplier to dispatch qualified maintenance personnel to repair.
Blood pressure cuff is not connected	Blood pressure cuff is not properly tied up or no blood pressure cuff	Tie up the blood pressure cuff
Blood pressure cuff air leakage	Blood pressure cuff, rubber pipe or connector is damaged	Check and replace the air leakage part, if necessary, inform the supplier to dispatch qualified maintenance personnel to repair.
Weak signal	Blood pressure cuff is too loose or the animal's pulse is too weak.	Use other method to measure the blood pressure.
Over pressure protection	The pressure is above the specified safety upper limit	Remeasure it again, if the failure continues, stop the NIBP measurement function, and inform the supplier to dispatch qualified maintenance personnel to repair.

Chapter 6 the Impedance Respiration Monitoring

This Monitor measures out the human body respiration movement waveform by the electrocardiogram lead line with the impedance method, this respiration movement waveform is used to calculate the human body respiration rate and analyses the human body respiration status.

The Monitor can measure three lead respiration waveforms according to the selection such as I, II, III, of which, the II lead is generally used. Because the different signal strength, the Monitor provides two kinds of gain adjustment functions such as automatic adjustment and manual adjustment, it can make the respiration waveform clear and without distortion by adopting appropriate gain.

6.1 The electrode connection

The impedance respiration method of this Monitor is to monitor the animal's respiration by the impedance variation between the two electrocardiogram electrodes. The respiration monitoring is realized via the electrocardiogram I, II or III lead and their corresponding electrodes, so no additional electrodes are required, and the installation of the electrode is very important. For some animal, because of their clinical situation, the horizontal expansion of the thoracic cage leads to the negativity thoracic cage internal pressure, under this situation, it is better to put the two respiration electrodes at the central line of the right side of the thoracic cage or the most active area on the left side of the thoracic cage when breathing, so that the optimum respiration wave is gained.

Here we suggest installing five lead electrodes for the respiration monitoring, as shown in Figure 6.1.

Notice

- ① **The respiration monitoring is not recommended to be used on the animal with large activity, because it will lead to mis-alarm.**
- ② **Symmetrically install white and red electrodes to get optimum respiration wave, avoid putting the liver and the ventricle on connecting line, which can avoid the false difference produced by the heart coverage or pulsation blood flow, this is very important for the neonatal.**

6.2 RESP monitoring inspection

- 1) Make skin preparation for the animal before installing the electrode;
- 2) Open the crocodile clips wide enough to grasp the skin.
- 3) Connect the power for the monitoring system.

Installing the electrode for the respiration monitoring system

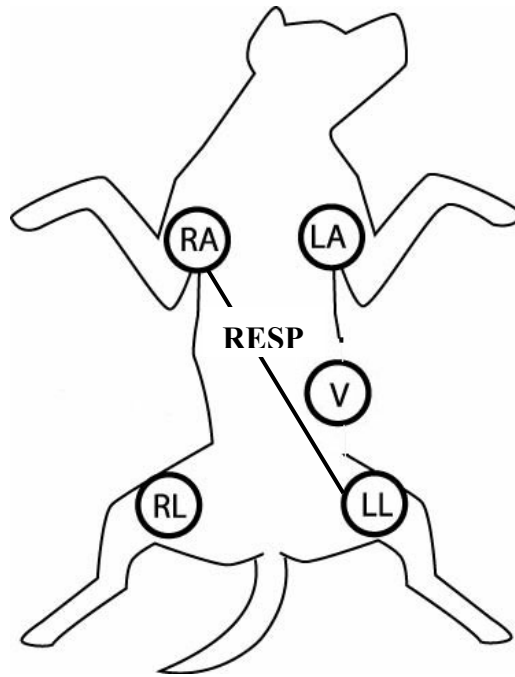



Figure 6.1

Notice

The white and red electrodes shall be installed on the cross to gain the optimum respiration wave. Avoid putting the liver and the ventricle on connecting line, which can avoid the false difference produced by the heart coverage or pulsation blood flow.


6.3 Monitoring interface

Impedance respiration monitoring interface and its menu are as shown in the Figure below:

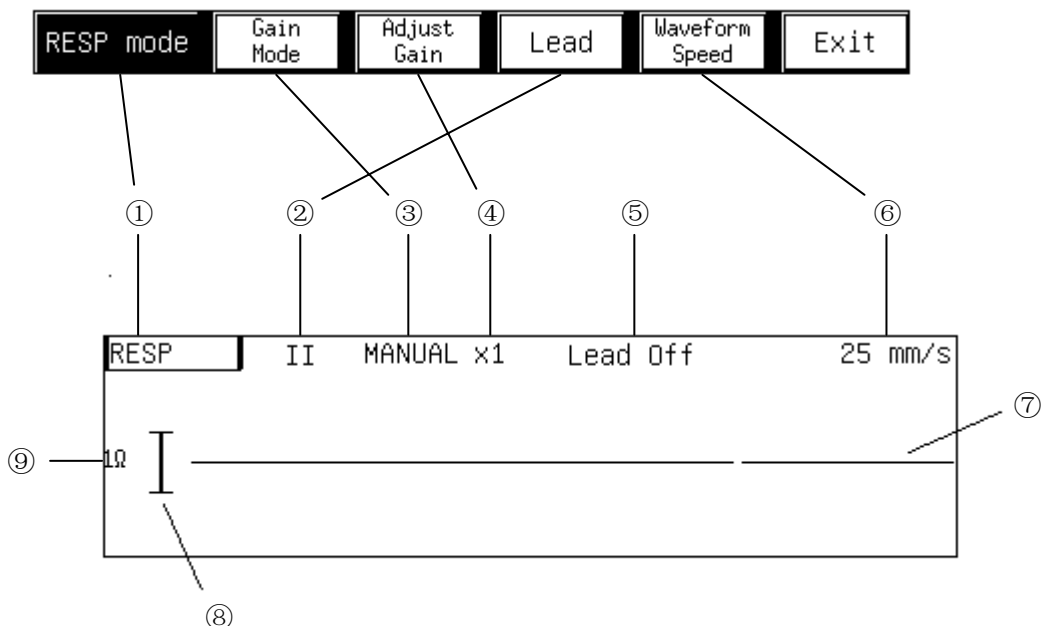


- ① Impedance respiration waveform label, select this label and the menu as the above Figure shows, the first knob "RESP mode" is used to switch the respiration measurement mode, which can be switched between the impedance respiration and nasal tube respiration.
- ② Lead type, its corresponding menu knob is the [Lead];
- ③ Gain mode, automatic or manual, their corresponding menu knob is [Gain Mode];
- ④ Gain, it has four grades such as x1/2,x1,x2,x4, their corresponding menu knob is [Adjust Gain];
- ⑤ Status prompting bar
- ⑥ Waveform speed is in four grades such as 6.25mm/s, 12.5mm/s, 25mm/s and 50mm/s, their corresponding menu knob is [Waveform Speed], at the same time, this knob adjusts the speed of the impedance respiration and nasal tube respiration.
- ⑦ Respiration waveform
- ⑧ Scale
- ⑨ The range corresponding to the scale. The unit is Ohm, the length and range of the scale is different from time to time along with different gains.

The respiration rate is displayed on the parameter area on the right side of the Monitor, as shown in the Figure below:

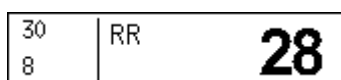

Select the respiration rate label, you can get following menu:

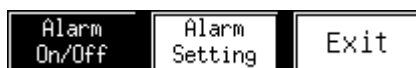

[Alarm On/Off]: allow and prohibit the respiration rate limit alarm.

[Alarm Setting]: check and set the respiration rate alarm setting.

6.4 Prompting messages

Lead off: it means that the current lead is not properly connected, when the lead line is not inserted on the Monitor or connected to the animal, this message shows.


Off line: it means the impedance respiration function is cut off internally, please check the settings in the menu key → [System Setup] → [Maintenance] → [Reserved 1][ECG Switch] → [I-RESP]. Please refer to the ECG on/off operation in the chapters and section about the ECG monitoring.

Running: except the abovementioned two conditions, the system will prompt "Running".

Chapter 7 Nasal Tube Respiration Monitoring

7.1 Monitoring interface

The nasal tube respiration monitoring interface and its menu are as shown in the Figure below:

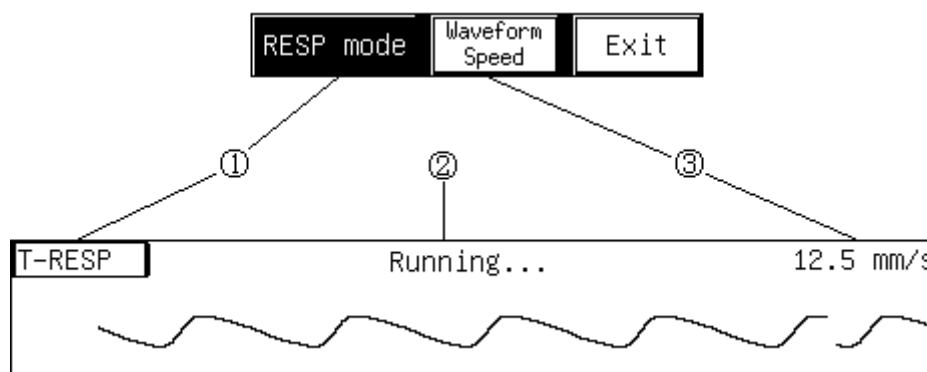

- ① Select this label and the menu as the above Figure shows, the first knob “RESP mode” is used to switch the respiration measurement mode, which can be switched between the impedance respiration and nasal tube respiration.
- ② Status prompting bar.
- ③ Waveform speed is in four grades such as 6.25mm/s, 12.5mm/s, 25mm/s and 50mm/s, their corresponding menu knob is [Waveform Speed], at the same time, this knob adjusts the speed of the impedance respiration and nasal tube respiration.

The respiration rate is displayed at the position of nasal tube respiration, they use the same position, please refer to the relevant parts concerning the impedance respiration monitoring.

7.2 The status prompting messages

Running: the system always shows “Running”.

Chapter 8 Temperature Monitoring

8.1 The installation of the temperature sensor

The G3D Multi-Parameter Animal Monitor can use two temperature sensors simultaneously, stick the body temperature sensors on the animal's body where ought to be measured, and two body temperature datum are measured out, and gain the body temperature difference.


Caution

- ⚠ Before the monitoring, check if the sensor cables are properly connected, pull the temperature sensor cable from the channel 1 jack, the error message "T1 sensor is off" is displayed on the screen and the alarm sound is produced, the other channels are similar to this.
- ⚠ Be cautious to deal with temperature sensors and cable. When it is not in use, the sensors and cable ought to be ringed loosely, if the cable cores inside them are too tight, which will lead to mechanical damage.
- ⚠ The temperature instrument calibration shall be carried out once a year (or carried out according to the time specified in the hospital procedure). If the calibration is needed, please contact the manufacture.

Notice

- ⚠ The disposable temperature sensor only can be used once.
- ⚠ During the monitoring process, the temperature measurement instrument will self-check every one hour, the self-check lasts for two seconds, it will not influence the normal operation of the temperature the Monitor.

The temperature measurement result is displayed on the parameter area of the Monitor, as shown in the Figure below:

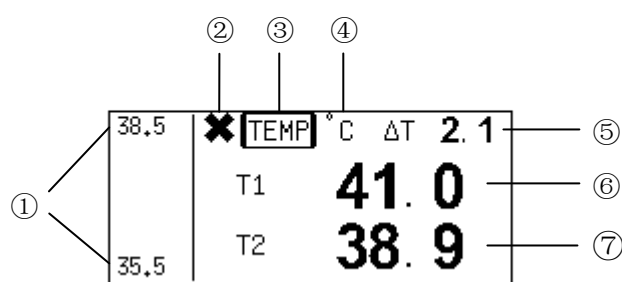

- ① Temperature alarm upper/lower limit.
- ② Temperature alarm off sign.
- ③ Temperature parameter label. Select this label, the body temperature parameter menu shows, which will be introduced later.
- ④ Temperature unit. It means the unit of the value.
- ⑤ Temperature difference. Temperature difference is the absolute value of the difference between the temperature1 and temperature2.
- ⑥ Temperature1, measurement range: 25.0-45.0°C.
- ⑦ Temperature2, measurement range: 25.0-45.0°C.

8.2 Body temperature parameter menu




【Alarm On/Off】 allow or prohibit temperature parameter alarm.

【Alarm Setting】 check and configure temperature alarm setting.

In this Monitor, tow temperature channel use one set of alarm limit and alarm on/off. Temperature difference has no alarm function.

【Mode】 selects display content, it has three combination forms, as shown below:

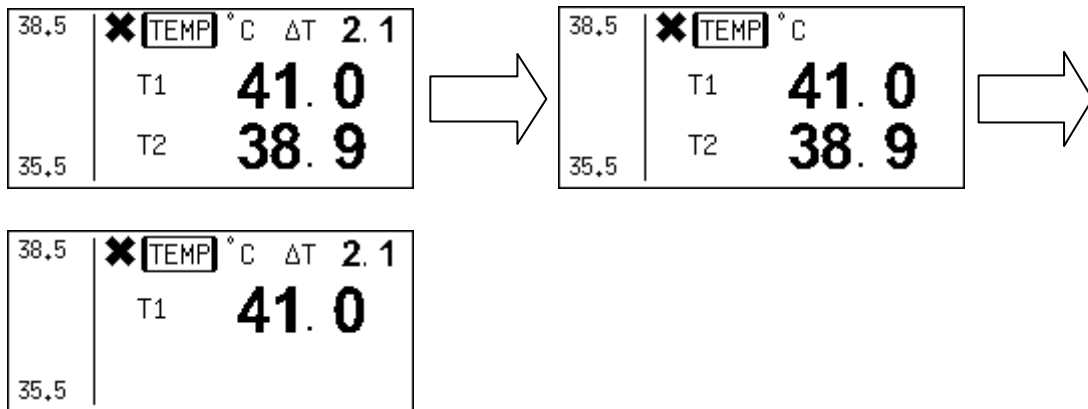

The doctor can select the display content as required.

8.3 Maintenance and cleaning

Caution

⚠ When washing the sensor connecting to the Monitor, turn off the power and cut off the AC, please.

Reusable temperature sensor

- 1) The heat for the temperature sensor can not be over 100°C(212°F). It can only endure the temperature of 80°C(176°F)~100°C (212°F).
- 2) The probe cannot be sterilized with the vapor.
- 3) It only can be sterilized with the detergent containing the alcohol.
- 4) When washing the sensor, hold the front end with one hand, the other hand clean the sensor toward the connector with a wet downless cloth.

Notice

⚠ If what you used is disposable temperature sensor, it cannot be repeatedly sterilized and reused.

⚠ In order to protect the environment, disposable temperature sensor shall be recycled or properly disposed.

Chapter 9 IBP Setting and Measuring (Optional)

9.1 Waveforms displaying


Two channels of waveforms are displayed in a same/different waveform area, differentiated by color. When the two curves overlapped, the color of curve will be the color of channel 2.

The height of the waveforms relates to the pressure. And the ratio between height and pressure is specified. Maximum value and Minimum value are listed in the following table:

Maximum value on screen	Minimum value on screen
6 mmHg	-1 mmHg
10 mmHg	-1 mmHg
18 mmHg	-2 mmHg
30 mmHg	-3 mmHg
60mmHg	-6mmHg
80mmHg	-8mmHg
100mmHg	-11mmHg
120mmHg	-13mmHg
180mmHg	-20mmHg
240mmHg	-26mmHg
300mmHg	-33mmHg

The monitor can automatically select a proper ratio according to the phase of the waveform. It also provides a function to adjust the ratio manually. Here we call it two gain modes: auto mode and manual mode.

At the left side of the waveform, there is a scaled axis. It is the pressure axis of channel 1. And at the right side of the waveform, there is also an axis for channel 2. The axis looks like:

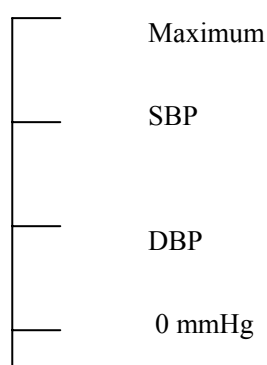

The bottom horizontal line indicates the location of zero pressure. If a curve dot is displayed higher than this line, it means the pressure of this dot is positive; accordingly, the pressure of the dot, which is lower than this line, is negative.

The top horizontal line indicates the maximum pressure of current display gain. The curve line where its pressure is higher than the maximum pressure will be clipped and becomes a section of horizontal line.

And the 2nd horizontal line counted from top is the location of SBP. The 3rd horizontal line is the location of

DBP. So via this scale ruler, operator can easily tell the pressure of the curve.

9.2 Module status display

The status of two channels is displayed at the line above the waveform. All the possible status messages are listed in the following table:

Message string	Meaning	Possible causation
No response	The monitor can not get data from this channel	When the monitor is restarting, this message indicates that the channel is not started yet, but the message should disappear in 5 seconds. The hardware of IBP module does not work properly. It needs to be fixed.
Zero process...	This message indicates zero operation is in progress.	This is a normal message. It will appear when zero is processing, and disappear immediately when the process is finished.
NOT zero	This message is to tell the operator that the channel did not perform zero processing	Operator has to perform zero process after the monitor is started. If not, the pressure measured is unusable.
(Auto/manual)	This is a normal working message. It tells the current gain mode and gain ratio.	The channel works properly.

9.3 Pressure value display

IBP measures SBP, DBP and MAP values. Pressure values of both channels are displayed at the right of the main screen. These values can vary in the range from -60mmHg to 300mmHg. Values not in this range will be displayed as invalid pressure, in the form of '---' or '--.-'.


But there are more reasons if these pressure value is displayed as invalid value, they are:

1. The actual pressure is not in the proper range from -60 mmHg to 300 mmHg.
2. The channel did not perform zero process.
3. The cable is disconnected or not connected firmly.
4. The channel is not calibrated properly.
5. Zero in processing.
6. Monitor is starting.

The pressure values are displayed in the same unit as the NIBP. So if operator wants the IBP pressure values are displayed in another unit, he has to use the menu of NIBP.

Menu operation guide

Let menu disappear, then use the rotating mouse to select a target IBP channel, press down the rotating mouse, we will get a menu which operates on the specified channel. In this menu, we can change the label of the channel, process zero, and select the display gain policy. Functions of these menu items are describing one by one below.

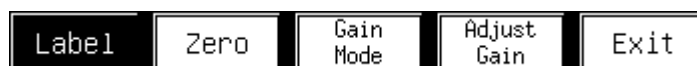

Select LABEL

Once the menu button is pressed, the label of the channel will be changed. We can see the changing of the label displaying near the pressure values and the one displaying in the status area. The label recycles in the order listed below:

Label	Customary meaning
ART	Artery pressure
CVP	Center vena pressure
RVP	Right ventricle pressure
LAP	Left atria pressure
RAP	Right atria pressure
PAP	Pulmonary artery pressure
ICP	Encephalic pressure
LVP	Left ventricle pressure

After the selection, the label will be saved in 5 seconds. The monitor will reload and assign the saved label to the corresponding channel at the next start.

These labels are provided as a sign to the channels only. It does not take any other effect other than displaying.

ZERO

In the first time we press this menu button, the corresponding channel will perform zero process. And press the menu button again will terminate the zero process. The system will treat the detected absolute pressure value as the relative zero point. This zero point affects both the pressure values and the waveforms of the channel.

Before performing a zero process, operator should connect the sensor and the pipes, and then plug the cable to the monitor. Make sure the other end of the pipe is open to the air and keep the pipe stay statically. Then it is ready for a zero process. The operator has to decide when the processing will finish. If the pipes are shaking, operator must wait until the pipes get static. If the pipes are static, the process can be done very fast.

Operator has to perform zero process after the monitor started. And before each measurement and after each time the pressure sensor or pipes is changed, zero process must be performed. It's recommended to perform zero process at least one time every day.

GAIN MODE

The menu button is to select the display gain mode for the channel. The gain mode switch between auto and manual when press this menu button. This function is disabled if the channel did not zero.

ADJUST GAIN

If the gain mode of the channel is manual, we can use this menu button to increase/decrease display gain of the channel. The clockwise for decrease, whereas for increase. This operation zooms in the waveform in vertical direction. This function is also disabled if the channel did not zero.

CALIBRATION

Calibrating operation is used to calibrate the sensibility of the pressure sensor. IBP module needs to change pressure sensor frequently. And there must be different sensibility between two sensors. So we need to calibrate when the pressure sensor is changed or the temperature changed. The result of calibration affects both waveform and pressure values.

Before calibrating, zero process must be done. The system will disable calibrate function if the channel did not zero.

Calibration can be done by the following steps: Connect the sensor pipe with a sphygmomanometer, then pump the sphygmomanometer until the azoth column rise to about 200mmHg. At the same time, the monitor will display the pressure it detected, we can see the pressure may not equal to the value of sphygmomanometer. If the azoth column is static, operator push the menu button→System Setup→Calibration→IBP Calibrate, enter the calibrate function of the channel, and input the correct value. In the calibrate function, it display the pressure value on the screen, then the operator can use the rotating mouse to change the value. Then press the rotating mouse to confirm the operation, or press the menu key to cancel. Once the new value is inputted, the system will calculate a new ratio and apply it immediately. The new ratio will be saved in 5 seconds, and will be load at the next start.

Chapter 10 Trend Analysis

The trend data is the animal's data collected along with the time accumulation, the trend analysis is animal status chart displayed according to the trend data; the trend graph or trend table can be prepared according to the parameters selected by the user. Except the noninvasive blood pressure, they all are consecutive curves, the trend graph (table) can provide time zoom and range zoom to meet the needs of observing trends map of the animal's status; at the same time, the corresponding trend data list is provided for detailed analysis.


10.1 Main interface trend analysis

Main interface trend graph (table)


Now begins to introduce the main interface trend graph (table): The main interface trend graph (table) displays the latest trend data distribution in-real time, and the parameters can be selected according to the user's intention. Turn to the time concerned and adjust the time zoom spans.

- Display the trend graph (table):

Press menu key and enters into Display, then click the Format Setup, in the display channel selection window select the trend graph (table), after determined, you can see the corresponding trend graph (table) in the main interface (as below).

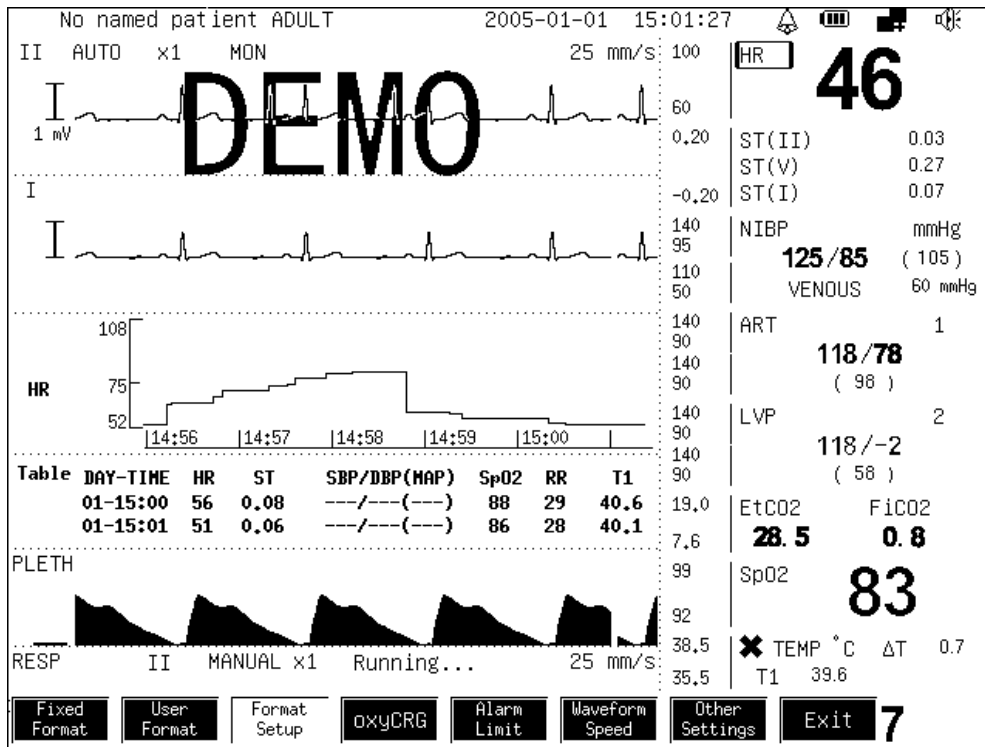

- Trend graph distribution introduction:

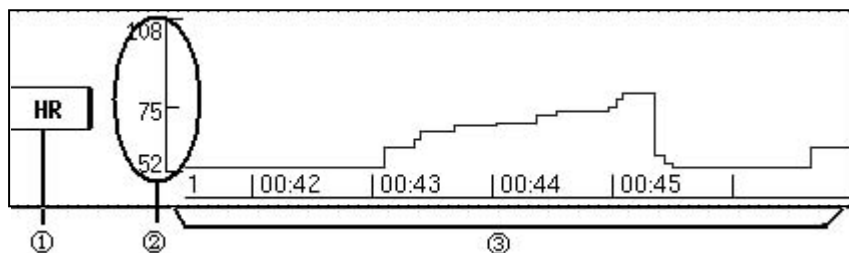

- ① **Parameter label:** The menu entrance for operation on the trend graph.
- ② **Measure scale:** Mark the value scale on the trend graph.
- ③ **Trend graph display area:** The trend graph is displayed on the upper part, the scale in the lower part is the time scale.


Trend graph display area: in the ordinary status, the trend graph of the latest time span is displayed in this area, the data is more than one screen, the graph automatically rolls left; you can use the “Last Page” to check the early stage trend data chart, use the “Next Page” to check the latest trend data graph, the color of the graph is determined by the parameters, the void values is brown.

Operational specification:


1. Label: click the parameter label and enter into the operational menu.



2. Parameter selection menu: the selection is as shown in the Figure below: any parameter; the corresponding trend graph is displayed in the trend graph display area.

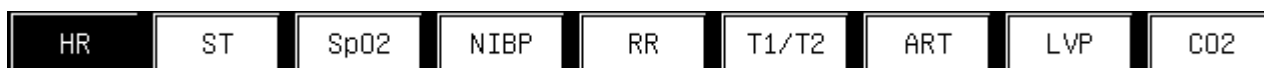

3. Adjust the time quantum of the trend graph: operate the [Last Page], [Next Page] in the menu to adjust the time, or turn the graph up or down.
4. [Step]: Operate the “Step” in the menu and enter the function as shown in the Figure below: select the step menu, the graph will be compressed according to selected step, the new graph in the display area of the trend graph line shall be added according to this step.

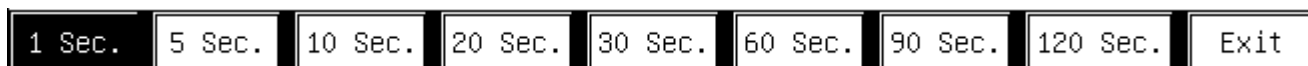

5. Range: The menu is shown as below, it is used to control the Y coordinate value range of the trend graph, the adjusted range will be stored in the Monitor, it will be applied at next start up; the range zoom has three adjustment methods such as alarm limit range, maximum range, manual adjustment, when the user set the zoom range, the trend data takes it as the upper limit/lower limit according to this range, the value over the limit is void value.



- Trend table distribution introduction

Table	DAY-TIME	HR	ST	SBP/DBP(MAP)	SpO2	RR	T1
	01-00:51	58	0.10	---/---(---)	89	30	40.8
	01-00:52	49	0.08	---/---(---)	85	28	39.9
	01-00:53	39	0.05	---/---(---)	80	25	38.9
	01-00:54	69	0.12	---/---(---)	94	32	41.9
	01-00:55	77	0.14	---/---(---)	98	34	42.7

①

②


- ① Menu operation: The menu entrance for the trend table operation.
- ② Data display area: In this area under the ordinary status, the trend data of the latest current time quantum is displayed in this area, when the data is more than one screen, it automatically rolls down or up, or use the [Last Page], [Line Up], [Line Down], [Next Page] to check the trend data of the time quantum concerned by the user.

Operational specification:

1. Enter into operation menu:




2. Adjust the time quantum of the data: adjust the time quantum of the time forward or backward data by using the knobs such as [Line Up], [Line Down], [Last Page], [Next Page];
3. Adjust the step: enter the step selection menu by operating the “Step” in the menu, which is shown in the Figure below: Select the step menu, the data in the trend table is compressed according to the selected step, the new data in the trend table will also be added according to this step.

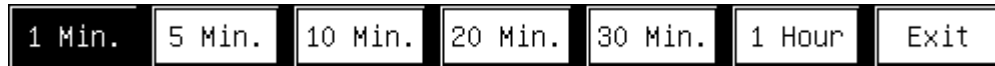

10.2 The trend analysis of the trend window

Now we introduce the popup trend window, it enters into the historic database from the system menu, and enters the trend window by clicking the trend display; the popup trend window can analyses the trend of one parameter or several parameters via the trend graph or trend table; it can conduct time zoom and manual adjustment value range according to the user’s intention.


- Trend graph distribution introduction:

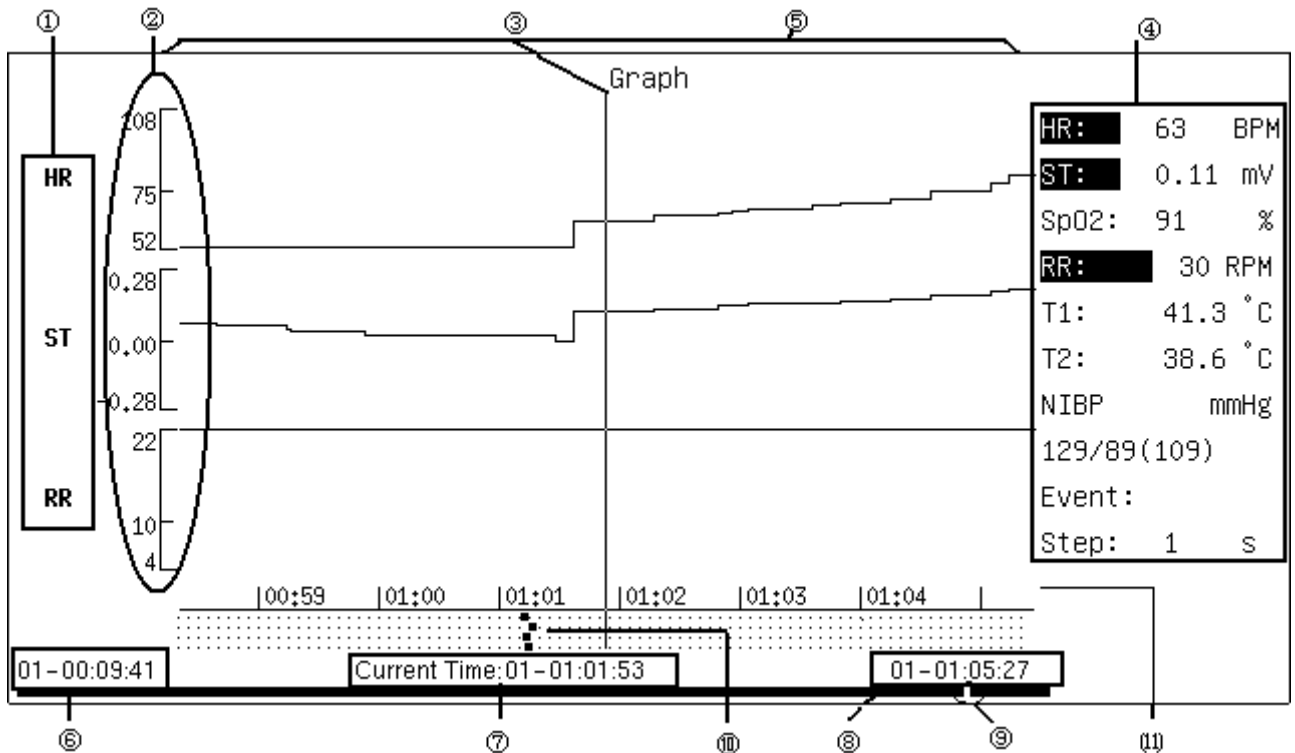


- ①**Parameter mark:** Mark the required drawing parameter selected by the user.
- ②**Scale:** Mark the value range of the trend graph.
- ③**Cursor:** Move this cursor by [Cursor].
- ④**Parameter table:** The corresponding value of the parameters at the current cursor time is displayed in this area, the inverse black text is the selected parameter, the part below the list is the area to display the event and step.
- ⑤**Trend graph display area:** Display the trend graph of the selected time quantum in the trend data, the color of the curve is determined by the parameter, of which, the void value is gray.
- ⑥**Starting time of the trend data:** It is used to mark the starting time of the trend data in the Monitor.
- ⑦**Cursor time:** It is used to mark the time value where the cursor is.
- ⑧**End time of the trend data:** It is used to mark the ending time of the trend data in the Monitor.
- ⑨**The mark of the cursor time value in the trend data:** It marks the cursor time value in the trend data with red point.
- ⑩**Event mark:** It is used to mark all event recorded in the current page.
- (11)**Time scale:** It is used to mark the time scale within the time range of current page.

Operational specification:


The operation menu is as following:

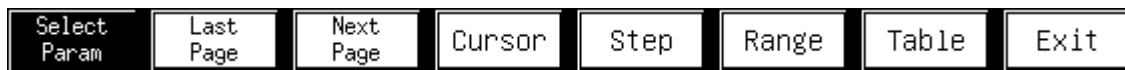

1. Select the parameter: The user select one to three parameter combination in the menu, the selected parameter is displayed with inverse positive figure in the parameter table, click the selected parameter again, it will be canceled, click the [OK], the corresponding trend graph will display in the trend graph display area.

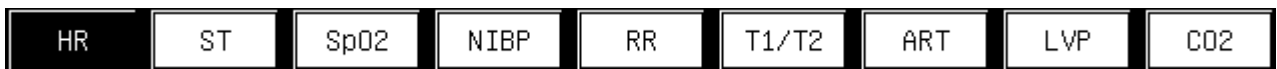

2. The graph time quantum: through the knobs such as [Last Page], [Next Page], the user can browse the trend graph forward or backward, until the required time quantum is founded.
3. [Cursor]: Click this knob and rotate, the cursor on the trend graph which move in accordance with the rotating direction, the time value and events where the cursor is will be displayed in the parameter table.
4. [Step]: select any time step as shown in the following Figure, the trend graph will compress the trend data according to this step, and display it in the trend graph display area.

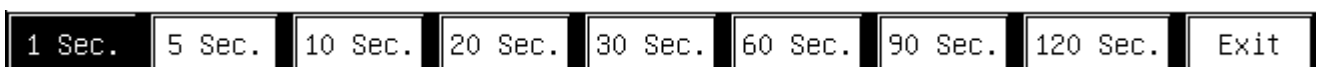

5. Range: The menu is as below, the value range used for control of the trend graph Y coordinate, the adjusted range will be saved in the Monitor, it is applicable in the next start up. There are three methods

for range zoom, including alarm limit range, maximum range, and manual adjustment, when the user set the time zoom range, the trend data will take this range as the upper limit/lower limit value, the value above this is voidable.



ALARMLIMIT: The alarm limit set by the system user is as the value range.

MAXRANGE: The default range of this type of Monitor serves as the value range.

MANUALADJUST: The user clicks and rotates the [MANUALADJUST] knob to adjust the upper/lower limit range, its value range is limited by the maximum range.


6. Trend table

Table			
DAY-TIME	HR	ST	SpO2
15-01:36	73	0.16	96
15-01:37	57	0.09	89
15-01:38	48	0.04	84
15-01:39	39	0.00	80
15-01:40	71	0.15	95
15-01:41	58	0.09	89
15-01:42	53	0.07	87
15-01:43	44	0.02	82
15-01:44	65	0.12	92
15-01:45	73	0.16	96

Line Up Line Down Last Page Next Page Step Graph Exit


In this trend table, the information in the trend graph is displayed in list, the parameter number and category are the same with the trend graph. In the trend table, the central time is the cursor time of the trend graph, in the trend table, the value of each parameter is date and time value (except NIBP), the NIBP is the first value in the current step.

Trend table operational specification:

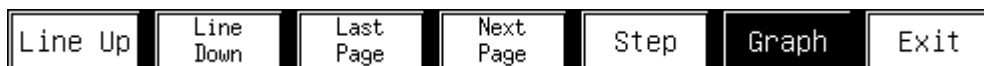

- Adjust the data display time range: Adjust the data display time range via the knobs such as [Line Up], [Line Down], [Last Page], [Next Page] to get the time quantum concerned by the user.

- Adjust the time step: Adjust the time span of each data line in the trend table by step.

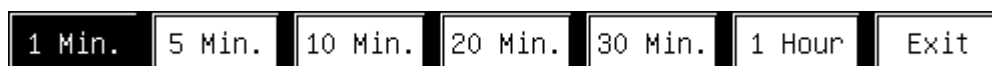

Chapter 11 the Waveform Freeze and Replay

Waveform freeze and replay, this function can freeze the current waveform of the Monitor, so that it is convenient for the doctor to check; turn several pages up or down, you can browse the current waveforms of several seconds, or save the required waveform for further analysis.

11.1 Waveform freeze operational specification

Press the freeze key in the control keyboard, the waveform in the screen is frozen, the following menu shows:



Browse the waveform: the user can turn the waveform forward or backward via the two knobs such as [Last Page] and [Next Page].

Save waveform: Press the [Save] knob to save the frozen waveform section for replay later, if the waveform is successfully saved, the following menu shows:




The waveform replay operational specification:

Enter the [Recall] menu via the [Review] knob, the first frozen waveform is displayed on the screen, the menu is as following:



Select frozen file: Press the [Select File] knob, the following interface shows:

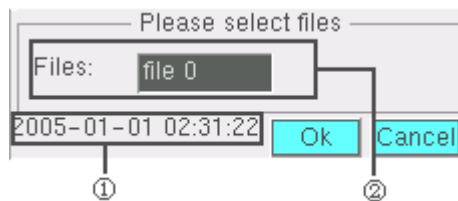

①The start time of the selected file;


②The selectable defrozen file.

Move the cursor to position ② you can select the corresponding time by rotating the knob, click [OK], the corresponding waveform is being drawn in the screen.

Browse waveform: The user can turn the waveforms forward or backward via the knobs [Next Page] and [Last Page].

Chapter 12 Event Record

Event record function, it is mainly used to record events for further reference, through the [Tools] menu, click the [Event], the following interface shows:

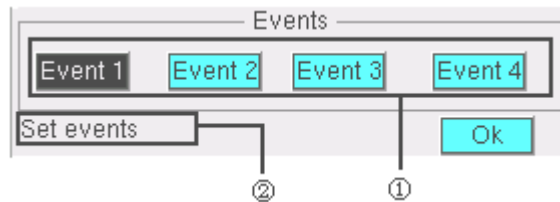

- ① It can be used to record the event numbers.
- ② Display that the event was successfully recorded.

Event record operational specification:

Record events: click the corresponding events recorded in the [Event1], [Event2], [Event3], [Event4], the position ② will display the event that was successfully recorded;

Browse events: you can see the event records at the corresponding time of the day in the event display area of the trend analysis.

Chapter 13 How to Use the Drug Calculator

The drug mixture used for intravenous transfusion and transfusion involves the information such as the drug dosage, infusion speed, dose, volume, and concentration. The drug calculator can calculate the unknown value via the known items to help you to control the drug infusion.

Caution Before using any drug, you shall check that whether you have selected the correct calculation unit or animal weight, if you have some problems, please consult the dispensary in your hospital.

It is always that the doctor in charge determines to use the drug and dosage. The drug calculator will calculate the dosage based on the input value, it cannot verify the validity of the calculated data.

1. Enter the drug calculator

In order to enter the drug calculator, press the menu key→Tools→Drug Calculator.

2. Unit

Drug unit

When the drug name is "Any Drug", you can select: g, mg, mcg, unit, k unit, m unit, mEq.

When you select the drug name, the drug unit is set by the drug calculator automatically, the user cannot modify it.

After the drug unit was set, the unit relating to the weight in the drug calculator changes correspondingly, the unit of the titration process table and titration table are keep pace with it.

Time unit

It has two selections such as min (minute), hr (hour). After setting the time unit, the units in the drug calculator changes correspondingly, the units in the drip table, titration table also keep pace with it.

3. Terms

- Amount

The total weight of the drug used by the animal within certain time.

- Volume

The volume of the mixed liquor formed by drug diluents and drugs.

- Dose/min, (dose/hr)

The drug quantity injected into the animal within unit time (per minute, hour)

- Dose/(kg*min), dose/(kg*hr)

The drug quantity per 1kg of animal weight injected in the animal within the unit time (per minute, hour)

$$\text{dose}/(\text{kg}*\text{min}) * \text{weight} = \text{dose}/\text{min}$$

$$\text{dose}/(\text{kg}*\text{hr}) * \text{weight} = \text{dose}/\text{hr}$$

- Rate

The volume of the mixed liquor injected into the animal body within unit time (per minute, hour), the unit is ml/min, ml/hr.

- Rate/kg

The volume of the mixed liquor per 1kg of animal weight injected into within unit item (per minute, hour), the unit is ml/min, ml/hr.

$$\text{Rate}/\text{kg} * \text{weight} = \text{rate}$$

- Time

The consumed time for drug transfusion (the unit is min, hr).

- Concentration

The concentration of the liquor form by the drug diluents and drug,

$$\text{concentration} = \text{amount} / \text{volume}$$

- Drip/min, drip/hr

The transfusion mixture within unit time.

- Drip/ml

The volume of each guttae dropped form the transfusion device, the unit is GTT.

4. Drug calculation execution

a) The operation of numeric input block

- The number is input according to the digit, when enter the numeric input block, the numeric input block selects the first digit at left, tune the rotating mouse the numeric input block selects each digit consecutively from left to right, when it reach the last, it jumps to the first digit repeatedly, press the rotating mouse and enter the digit selection status.
- The selection range for each digit is 0~9.
- When the digit is more than the display value, it displays ---.--, when it less than the display value, it displays 0.00.

- When it cannot display all digits, the round is adopted.

b) Calculation rules

- The relational expression of the drug weight

$$\text{amount} = \text{dose/min} * \text{time}$$

$$\text{dose/min} = \text{dose}/(\text{kg} * \text{min}) * \text{weight}$$

- The relational expression of liquid volume

$$\text{volume} = \text{rate} * \text{time}$$

$$\text{rate} = \text{rate/kg} * \text{weight}$$

- Concentration = amount / volume

c) Known item and calculation results

- Amount, dose/min, dose/(kg*min), volume, rate, rate/kg, time and concentration can be input as the known item or be output as the calculation result.
- The user input the known items (at least three) according to the calculation requirements, according to the input items; the drug calculator can automatically calculate the item which has computation relationship with it.
- The known item input by the user is expressed with digit in blue background, the calculation result is expressed with digit in gray background.
- When enters this window, all items are displayed as 0.00, all items can be input with value, when the user inputs the known item, the drug calculator can calculate and the result displays in-real time, then the calculation result is locked, the user only can modify the digital value in the input known item, modify the input item as 0.00 means cancel this input.
- The reset knob restores to the initial status of the input.

The screenshot shows a 'Drug Calculator' window with the following fields and values:

Drug	Any Drug	Weight	50	kg
Drug Unit	mg	Time Unit	min	
Drug Weight		Volume		
Amount	1 0 . 0 0	mg	Volume	1 0 0 . 0
Dose/min	1 . 0 0	mg	Rate	1 0 . 0 0
Dose/(kg*min)	0 . 0 2 0	mg	Rate/kg	0 . 2 0 0
Time	1 0 . 0 0	min	Concentr	0 . 1 0 0
Drip/min	1 5 0	GTT	Drip/ml	1 5
Reset		Drip Table		Titration Table
				Ok

d) Conduct calculation for any drug

- When the drug selection is any drug, the drug calculator only provides the calculation function for the transfusion liquid, while have no application range prompt for the dose, concentration.
- Input three items such as weight, drug unit, and time unit.

- According to the known situation, input the drug use into corresponding items, the drug calculator will calculate relevant item for you.

e) Conduct calculation for exact drug

- The drug calculator has been preset more than ten general drugs calculation, including Aminophylline and Amrinone Lactate (refer to Drug dosage range limitation table), when select the drug name, the drug's items such as concentration, amount, dose/min have corresponding limitation ranges, when the input item or calculation result is out of range, the drug calculator is expressed in pink digit.
- The drug dosage over range prompt for the drug calculator, it only prompts the user to pay attention to use appropriate dose of the current drug, the final transfusion dose and transfusion process is determined by the doctor.

Drug Calculator

Drug	Aminophylline	Weight	50	kg
Drug Unit	mg	Time Unit	min	
Drug Weight		Volume		
Amount	100.00	mg	Volume	100.0
Dose/min	1.00	mg	Rate	10.00
Dose/(kg*min)	0.020	mg	Rate/kg	0.200
Time	10.00	min	Concentr	0.100
Drip/min	150	GTT	Drip/ml	15

Drug dosage range limitation table

English name	Concentration	Speed	Dosage
Aminophylline	0.5~1mg/ml	<25mg/min	250~500mg
Amrinone Lactate	1~3mg/ml	5~10ug/kg/min	5~10mg/kg
Bretylium Tosilate	10mg/ml	1~50mg/min	5~10mg
Dobutamine	0.5~5mg/ml	2.5~40ug/kg/min	250mg
Dopamine	<3.2mg/ml	1~30ug/kg/min	10~20mg
Epinephrine	<64ug/ml	1~4ug/min	0.025~1mg
Heparin Sodium	20~40U/ml	15~20U/kg/h±10%	6000~20000U
Isuprel	2~4ug/ml	2~20ug/min	0.5~1mg
Lidocaine	<8mg/ml	1~4mg/min	1~2mg/kg
Morphine hydrochloride	0.5~5mg/ml	<2mg/min	5~15mg
Nitroprusside	0.1~1mg/ml	0.5~10ug/kg/min	50mg
Nitroglycerin	50~400ug/ml	5~200ug/min	5~10mg
Oxytocin	0.01U/ml	0.001~0.04U/min	2.5~5U
Procainamide	2~4mg/ml	1~6mg/min	10~15mg/kg

5. Display the drip table

The “Drip Table” shows how much liquid is transfused into the animal’s body at a glance, and how much time remains.

- The drip table displays the residual quantities of the drug and liquor at each time quantum when the user inputs the data into the drug calculator window.
- At the main window of the drug calculator, select “Drip Table” knob and enter this function.
- The left side of the widow shows the data (and calculation result) inputed by the user at the drug calculator window, the right side is equally divided into 15 parts according to the whole transfusion process, lists the amount and volume at each time quantum.
- All items in this window cannot be modified.
- In the table, the unit of the amount is the same as the input unit at the drug calculator main window, the volume unit is ml.
- Press the “OK” knob and exit the “Drip Table” window, return to the drug calculator main widow.

Aminophylline			Drip Table		
Weight	50.0	kg	Amount	Volume	Time
Amount	10.00	mg	0.67	6.67	0:00:39
Dose/min	1.00	mg	1.33	13.33	0:01:19
Dose/(kg*min)	0.020	mg	2.00	20.00	0:02:00
Volume	100.0	ml	2.67	26.67	0:02:39
Rate	10.00	ml/min	3.33	33.33	0:03:19
Rate/kg	0.200	ml/min	4.00	40.00	0:04:00
Time	10.00	min	4.67	46.67	0:04:40
Concentr	0.100	mg/ml	5.33	53.33	0:05:19
			6.00	60.00	0:06:00
			6.67	66.67	0:06:40
			7.33	73.33	0:07:19
			8.00	80.00	0:08:00
			8.67	86.67	0:08:40
			9.33	93.33	0:09:20

Ok

6. Display the titration table

At a glance, you can see the drug dosage accepted by your animal at different rates, the higher the rate is, the larger the steps are between the items of the table.

- The titration table shows that at the same liquid concentration, the dose/min has a coincidence relation with the rate.
- Select the “TitrationTable” knob at the drug calculator main window and enter this function.
- The left side of the widow shows the data (and calculation result) input by the user at the drug calculator window, the right side is equally divided into 30 parts according to the range in the benchmark item, it shows the coincidence relation between the dose/min and rate.
- In the list the user has two items such as ”Dose/min”(if the time unit selected in the drug calculator main window is hr, then this item and titration table list head both are displayed as” Dose/hr”) and “Rate”.

- Dose/min (Dose/hr), its range is the one to two items as the dosage input by the user in the drug calculator, calculate the titration table.
- Press the “OK” knob and exit the” Titration Table” window, return to the drug calculator window.

Aminophylline			Titration Table			
Weight	50.0	kg	Rate	Dose	Rate	Dose
Amount	10.00	mg	0.40	0.04	6.40	0.64
Dose/min	1.00	mg	0.80	0.08	6.80	0.68
Dose/(kg*min)	0.020	mg	1.20	0.12	7.20	0.72
Volume	100.0	ml	1.60	0.16	7.60	0.76
Rate	10.00	ml/min	2.00	0.20	8.00	0.80
Rate/kg	0.200	ml/min	2.40	0.24	8.40	0.84
Time	10.00	min	2.80	0.28	8.80	0.88
Concentr	0.100	mg/ml	3.20	0.32	9.20	0.92
			3.60	0.36	9.60	0.96
			4.00	0.40	10.00	1.00
			4.40	0.44	10.40	1.04
			4.80	0.48	10.80	1.08
			5.20	0.52	11.20	1.12
			5.60	0.56	11.60	1.16

Ok

7. Reset

Exit the drug calculator window or turn off the Monitor, the drug calculator data input by the user is still saved, when the user enters the drug calculator window again, the latest input shows. If the user wants to start a new calculation, please press the” Reset” knob in the drug calculator and clear the input data, the new calculation starts.

Chapter 14 The Printing Function of the Monitor (Optional)

For the Monitor, the printing function is optional, you determine whether to have the printing function according to the different configuration of the Monitor, which can print one real time waveform of the main lead ECG (main lead definition please refer to other Chapters and sections) or two real time waveforms of the main lead ECG blood volume, with two printing speed of 25mm/s, 12.5 mm/s; and provides timing printing and alarm printing.

14.1 The connection mode of the printer


The printer used for this Monitor is a special thermal printer, which adopts embedded installation mode, generally it is installed on the side of the Monitor. Do not install any other types of printer to this Monitor, or it will damage the printer. If you have any doubt, please contact the manufacturer immediately.

- When the printer is properly installed, turn on the Monitor, the indicator of the printer power is green, if the red error alarm indicator lights, it means that the printer is in paper lack status, do not print at this time, or it will lead to the damage of the printer.
- The printing paper is single side thermal printing paper, lightly press the printer knob, the cover opens, the side of the printing paper with temperature sensitive coating is smooth, this side shall be installed upside, close the printer door, the red error alarm indicator is off.
- Check the printer whether it's properly connected, if not, the printer setting popup block will show the information "Printer is not probed correctly!".

14.2 The configuration of the monitor

When setting the printer, please select the "Recorder Setup" menu in the main menu, popup the dialog box "Set Printer Parameters".

The printer setting interface is as follows:

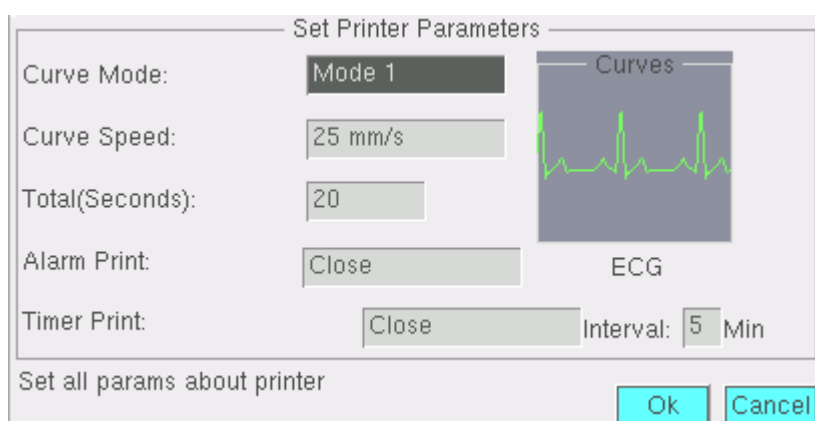

Curve mode is used to select the curve channel numbers and content which need to print, in "Mode 1" it prints one electrocardiogram, namely: the main lead ECG; in "Mode 2" it prints two waveforms: the main lead ECG+blood volume waveform.

The printer paper speed can be regulated via the "Curve Speed", it has two grades to be selected: 25mm/s,

12.5mm/s.

“Total” is the total printing length when printing waveform by pressing print key. The unit is second, the minimum length is 5 seconds, the maximum is thirty seconds.

“Alarm Print” is to print all the parameters when a parameter alarms, including heart rate, SpO₂, blood pressure, respiration rate and body temperature.

“Timer Print” is only to print the parameter. The default settings for the “Alarm Print” and “Timer Print” are off.

- **Notice:** when the printer is in working status, it cannot be configured.
- **Notice:** When the timing time arrives, if the printer is in working status, the Monitor will give up this parameter printing.

14.3 Waveform and parameter printing

Press the print key and start up the waveform printing, when reaching at the setting print time, it stops waveform printing and start to print the parameter, if the printer does not stop printing, press the print key and stop the current printing in advance.

- **Notice:** when the thermal recorder start up to printing the waveform, you can not forcefully open the printer’s cover when the recorder is still in work, you can take the printing paper out unless you stop the recorder.

If the printing paper in the recorder is used up, you must install the printing paper again, then restart the recorder, you could not start the printing function without the printing paper.

Chapter 15 CO2 Setting and Measuring (Optional)

15.1 Introduction

CO2 monitor is to monitor the respiration of animal by detecting the concentration of CO2 generated during respiration. The maximum concentration of CO2 at the end of exhalation is called *End-Tide CO2* (ETCO2). The minimum concentration of CO2 at the end of inspiration is called *Inspiration CO2* (inCO2). CO2 is generated by cells in the body during metabolizing, and is breathed out via breath system. The concentration of CO2 breathed out from lung reflects directly the situation of metabolizing and breathe system. If the concentration of CO2 is high, it means that metabolizing is excited, such as blood poisoning or acute fever. If the concentration of CO2 is low, it is commonly because the output ability of heart is weak, or the heart stop beating, or the lung soak for the flux of blood is low or oxygen it carried is not enough. Monitoring CO2 is used to warn the doctor of the abnormal of breathe and metabolizing for hocused animal.


The concentration of CO2 is represented as a pressure level, with 'torr' or % as its unit. Generally, the acceptable value is 38mmHg (5%) when air pressure is 760mmHg. The concentration of CO2 varies very fast from 0% to 5% normally. To detect the concentration of CO2 accurately, the monitor has to be very sensitive.

15.2 CO2 measuring

Firstly plug the air-water separator to the socket. Then push the bottom of the separator carefully to its position. Then connect one end of the breath pipe to the separator and the other end to the animal.

Secondly, we use the menu of monitor to enable the CO2 measure function as the follows:

Select **menu key** -> **Display** -> **Format Setup**

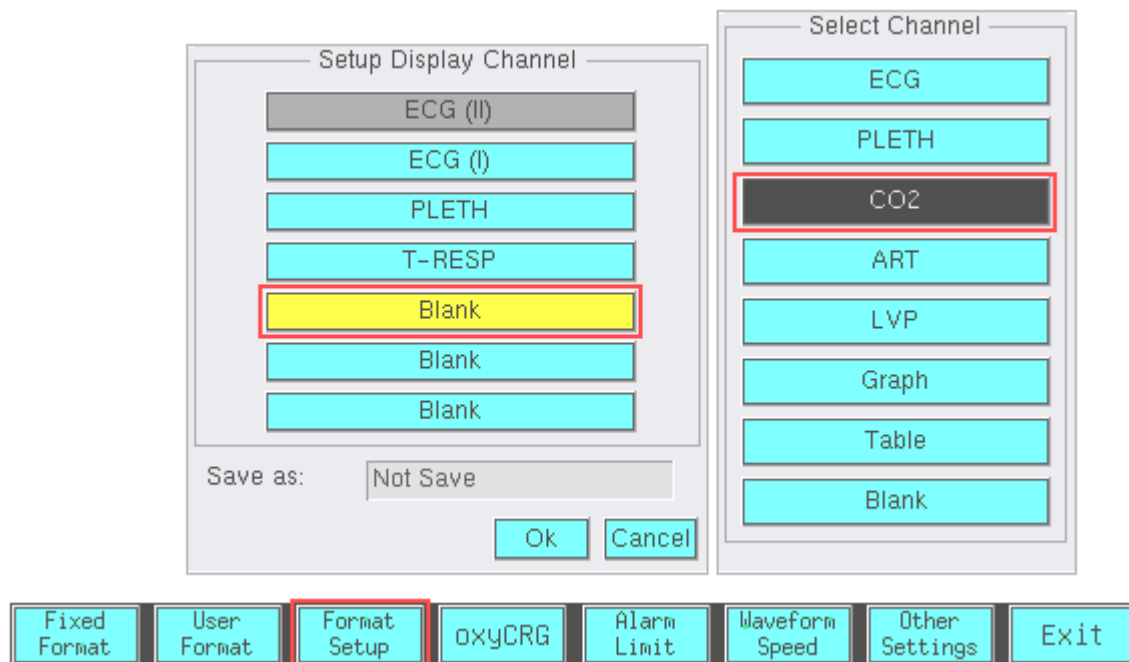


When the CO2 channel displays on the screen, the Monitor works in the CO2 respiration monitoring status. At this time, the value of ETCO2 and inCO2 are displayed on the parameter displaying area. (Notice: if CO2 channel isn't displaying, the CO2 respiration monitoring is closed.)

Important Note: please close CO2 channel (means do not display CO2 waveform on screen) when disconnect

CO2 accessories with animal and do not use CO2 function. If not, the motor in the CO2 module keeps running. Doing this can prolong the life of CO2 module.


15.3 Monitoring interface

CO2 monitoring interface and its menu are as shown in the Figure below:

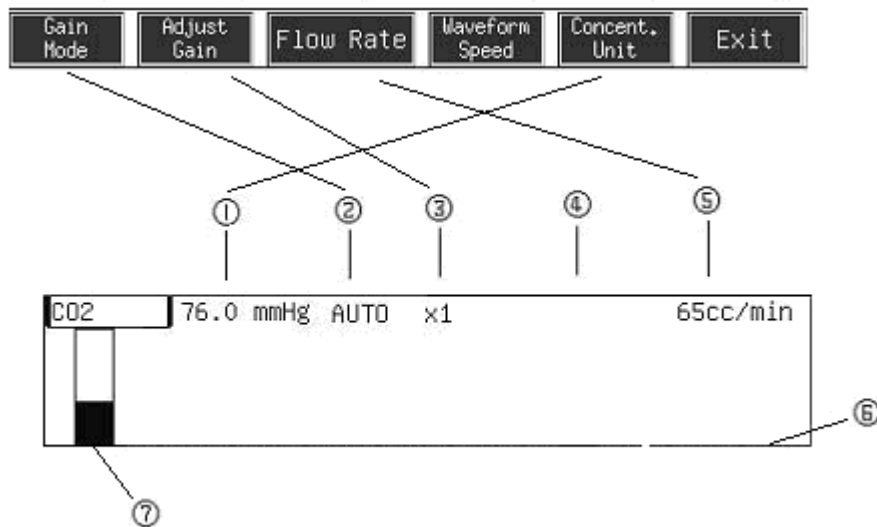

- ① Concent. unit, it has three grades such as 76.0mmHg, 10.13Kpa, 10.00%, their corresponding menu knob is [Concent. Unit].
- ② Gain mode, automatic or manual, their corresponding menu knob is [Gain Mode].
- ③ Gain, it has five grades such as x1/2,x1,x2,x4,x8, their corresponding menu knob is [Adjust Gain].
- ④ Status prompting bar.
- ⑤ Flow rate,it has three grades such as 65cc/min, 100cc/min, 150cc/min. When monitoring CO2, doctor can adjust the detect flux basing on the situation of the animal. Their corresponding menu knob is [Flow Rate].
- ⑥ Respiration waveform.
- ⑦ Scale.

The EtCO2 and inCO2 are displayed on the parameter area on the right side of the Monitor, as shown in the Figure below:

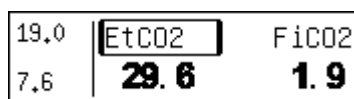

Select the EtCO2 label, you can get following menu:





[Alarm On/Off]: allow and prohibit the EtCO2 and inCO2 limit alarm.

[Alarm Settings]: check and set the EtCO2 and inCO2 alarm setting.




[Concent. Unit]: concent unit setting.

Chapter 16 Packaging and Label

16.1 Packaging frontispiece

		ANIMAL MONITOR
Model: G3 Quantity: 1pcs Date:		
Power Range: ~90V-265v, 50/60Hz, 80VA		
Manufacturer: General Meditech, Inc.		
Tel.: 86-755-26546289		Post Code: 518057
Website: http://www.szmedtech.com		
Address: South Office 4/F, Kezhi 1st Rd. West, Science Park Nanshan, shenzhen, Guangdong P.R. China		

16.2 Packaging profile

			
HANDLE WITH CARE			
Net Weight: _____ kgs			
Gross Weight: _____ kgs			
Dimension: 425X300X350mm			
C/No.: _____ of			

16.3 Product nameplate

ANIMAL MONITOR	
Model: G3	S/N: <input type="text"/>
Standard: YZB/粤 0320-2003 Q/SZGR1-2004	
Registered No.: 粤食药管械(准)字 2004 第 2210391 号	
Power Range: ~90V-265V, 50/60Hz, 80VA	
General Meditech, Inc.	